

Cadenza CW-C18

Unison UK-C8



Unison UK-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Carbazochrome ( Hemostatic )

カルバゾクロム ( 止血薬 )

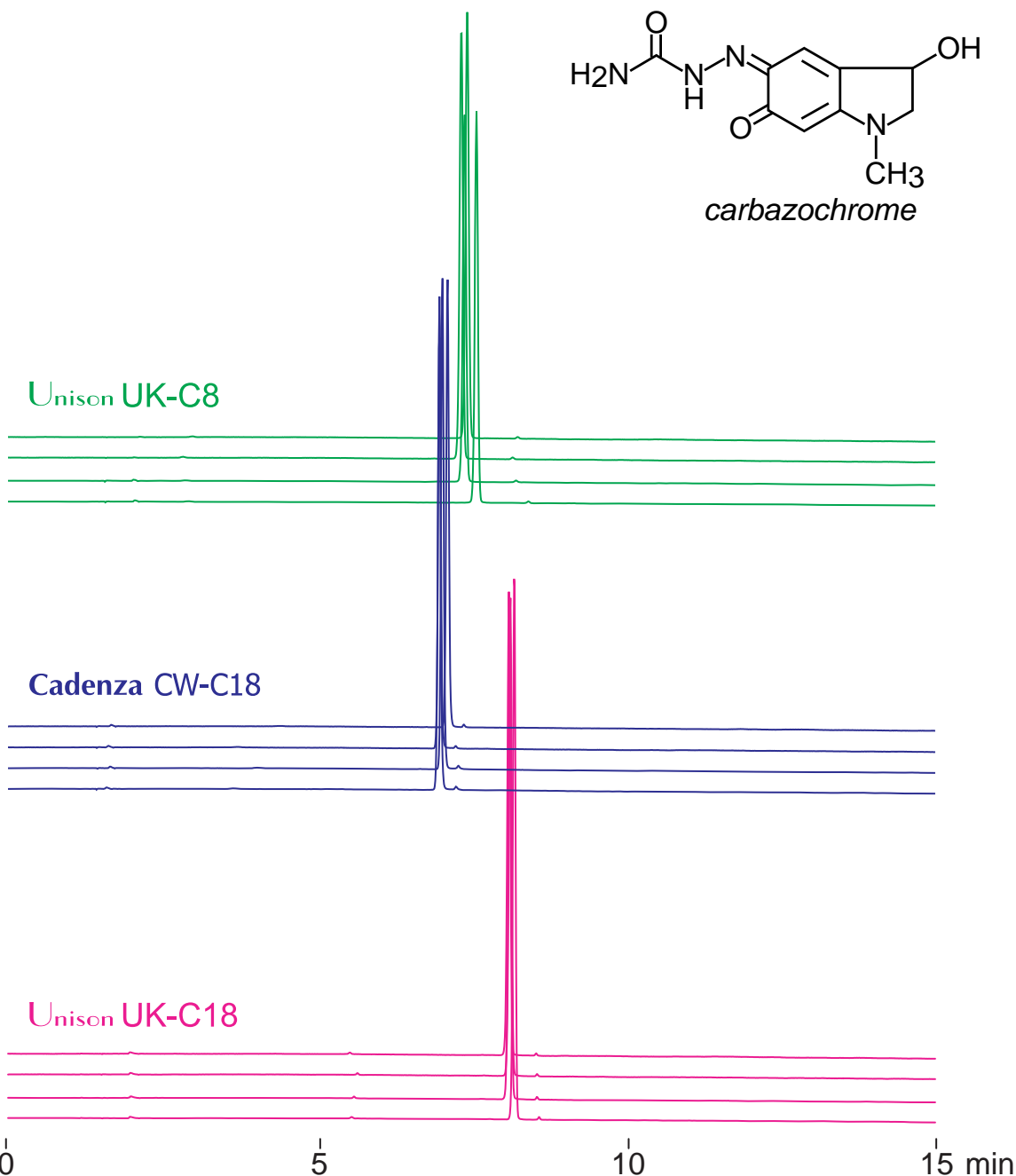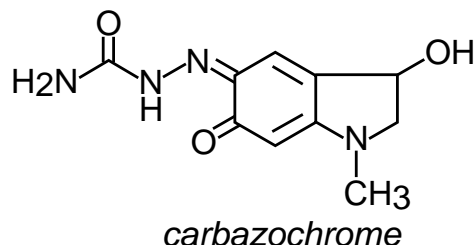

100 x 4.6 mm

A: 100mM ammonium acetate

B: 100mM ammonium acetate / acetonitrile = 10 / 90

0-50%B (0-15min)

1.0 mL/min (12MPa), 37°C, 360nm

10µL (0.05mg/mL, 100mM ammonium acetate)