

Cadenza CW-C18

Unison UK-C8


Unison UK-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Carbazochrome (Hemostatic)

カルバゾクロム (止血薬)

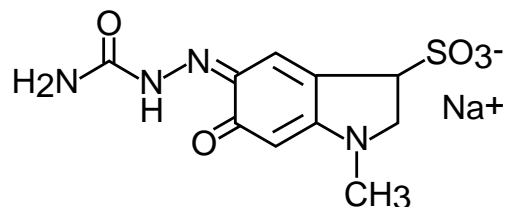

carbazochrome sodium sulfonate

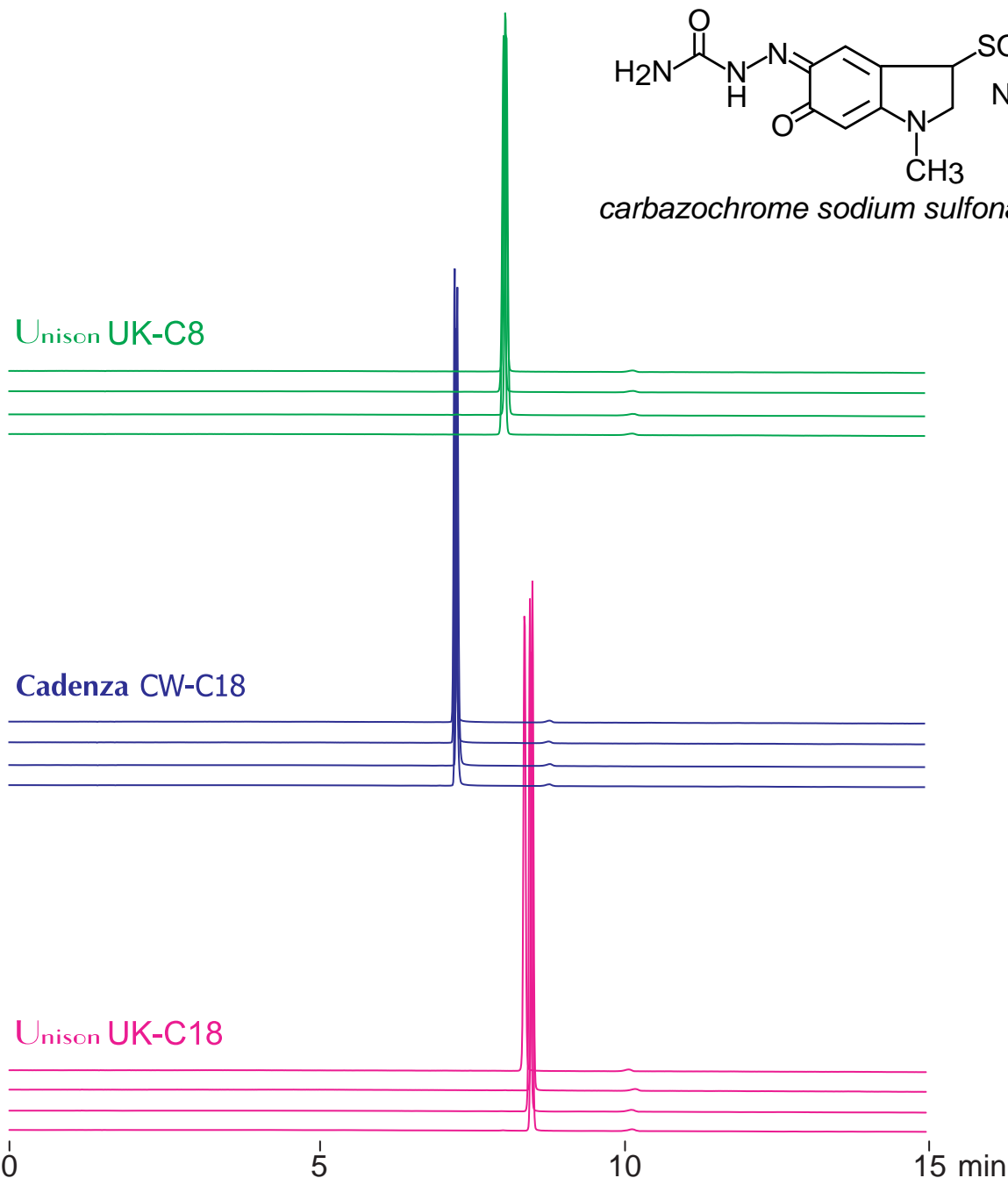

100 x 4.6 mm

A: 100mM ammonium acetate

B: 100mM ammonium acetate / acetonitrile = 10 / 90

0-50%B (0-15min)

1.0 mL/min (12MPa), 37°C, 360nm

10 μ L (0.05mg/mL, 100mM ammonium acetate)