



Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Ifenprodil ( Vasodilator )  
 イフェンプロジル ( 末梢血管拡張薬 )

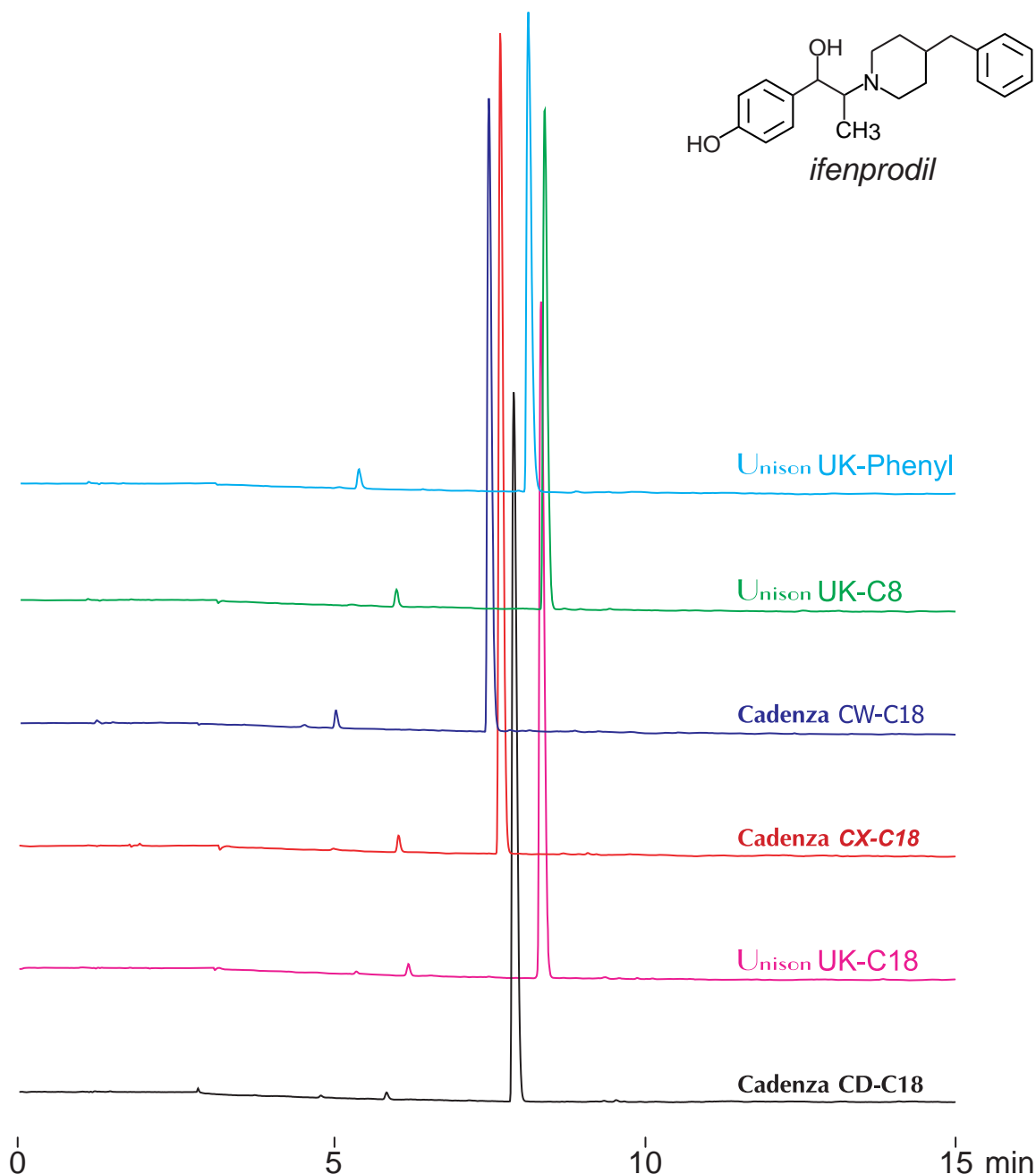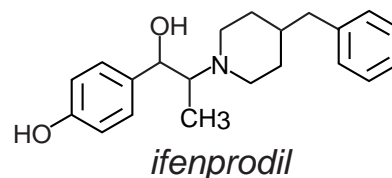

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-100%B (0-15min)

1.0mL/min (8MPa), 37°C, 280nm

10µL (1.0mg/mL, 100mM ammonium formate)