


Unison UK-Amino

250 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Piperidone related compounds

ピペリドン関連化合物の LC-MS 分析

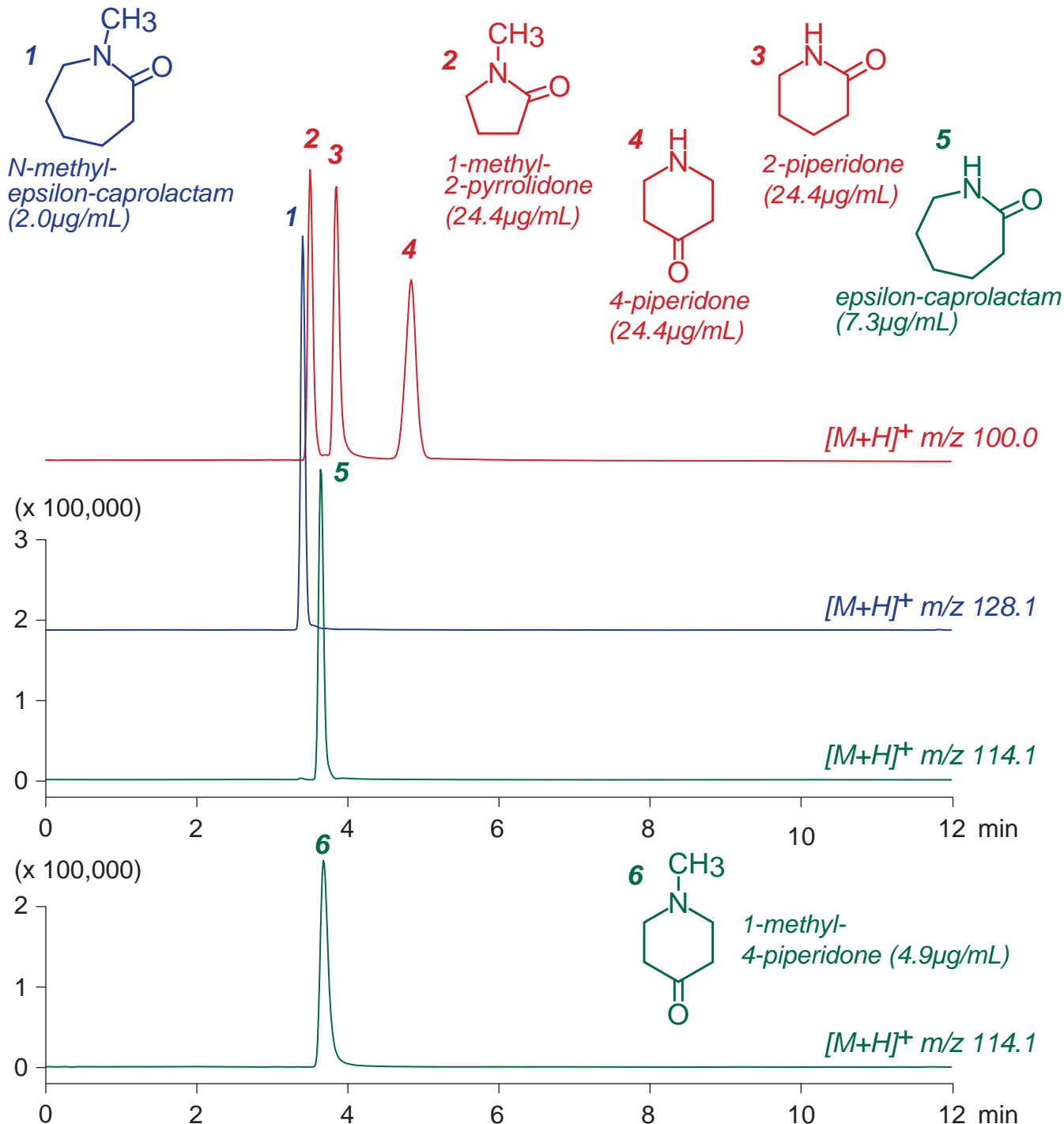

Unison UK-Amino, 250 x 3 mm

A: acetonitrile

B: 100mM ammonium formate

5-100%B (0-10min), 100%B (10-12min)

0.4mL/min (5MPa), 60°C, Single Quad. MS (ESI)

2 μL (2.0-24.4 μg/mL, ACN / 100mM HCOONH₄ = 50 / 50)