
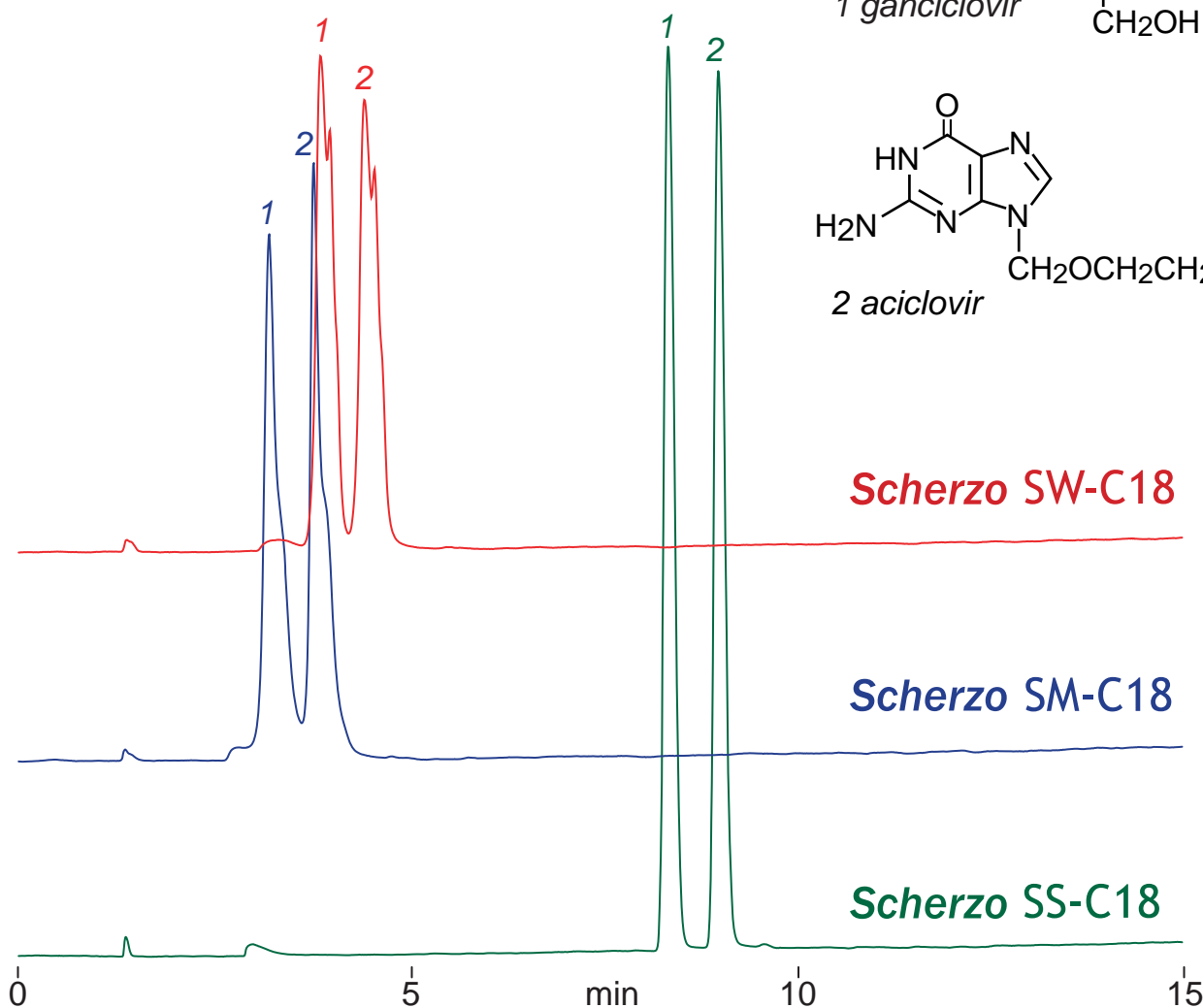
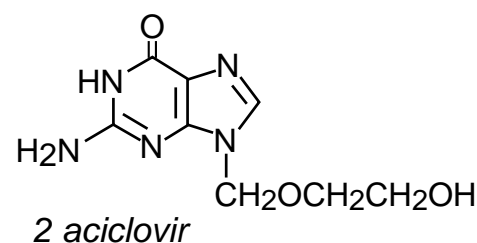


Scherzo SW-C18
 Scherzo SM-C18
 Scherzo SS-C18

150 x 4.6 mm

HPLC COLUMN APPLICATION

Ganciclovir, Aciclovir (Antiviral)
 ガンシクロビル, アシクロビル (抗ウイルス薬)

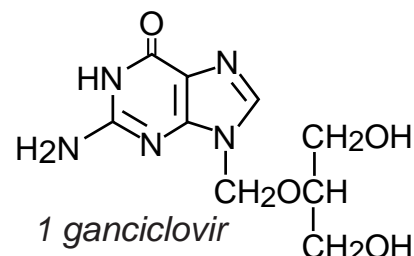

100 x 4.6mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 50 / 50 / 0.5

0-70%B (0-15min)

1.0mL/min (13MPa), 37°C, 260nm

10µL (0.1mg/mL, water / formic acid = 100 / 1)