

Unison UK-Phenyl Unison UK-C18



Unison UK-C8

Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Aminophenazone (Anti-inflammatory)

アミノピリン (ピラゾロン誘導体, 消炎鎮痛薬)

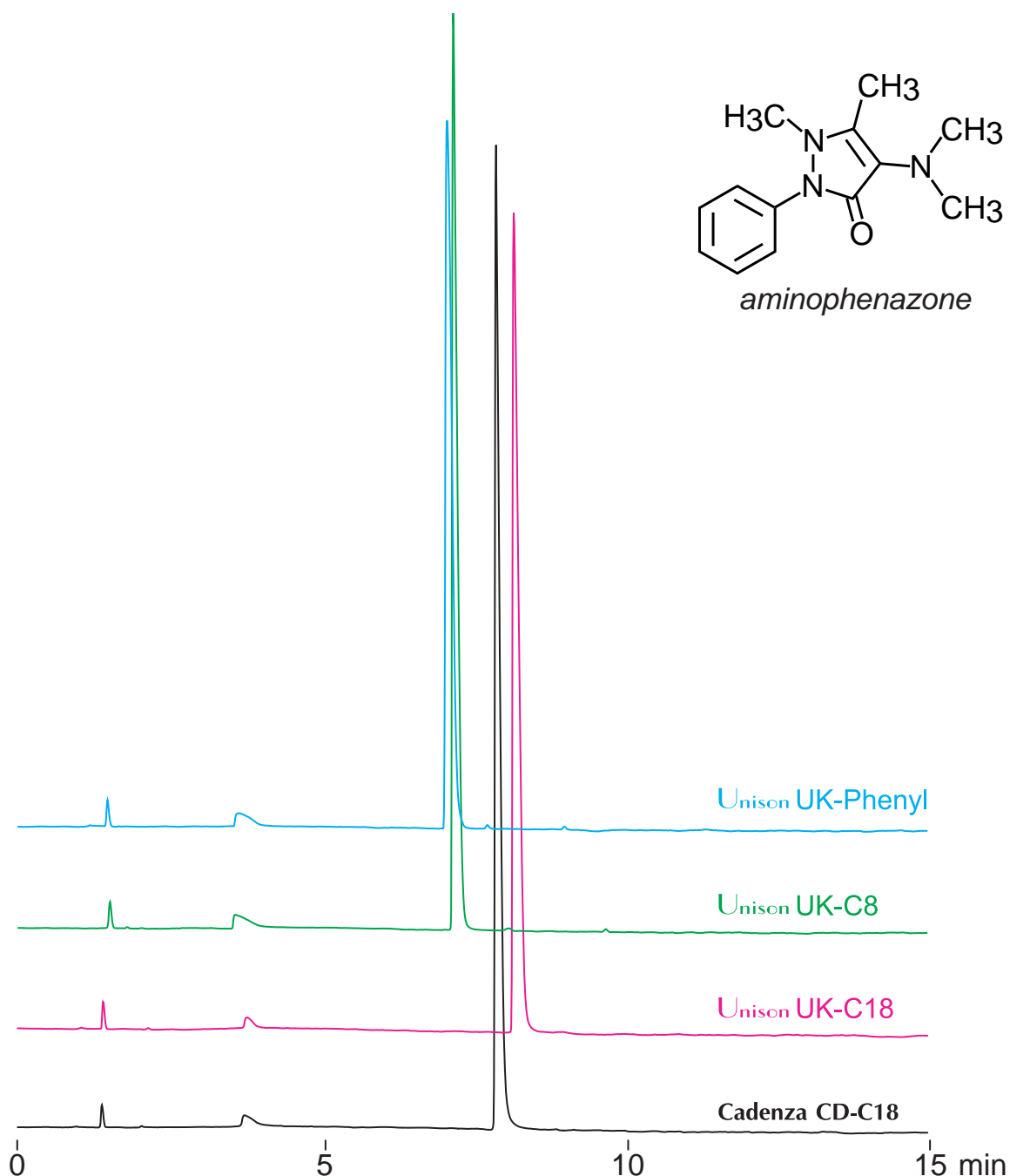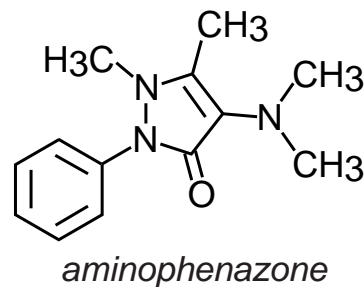

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-50%B (0-15min)

1.0mL/min (12MPa), 37°C, 260nm

10μL (0.1mg/mL, water / formic acid = 100 / 1)