

Unison UK-Phenyl Cadenza CW-C18



Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Aminophenazone (Anti-inflammatory)

アミノピリン (ピラゾロン誘導体, 消炎鎮痛薬)

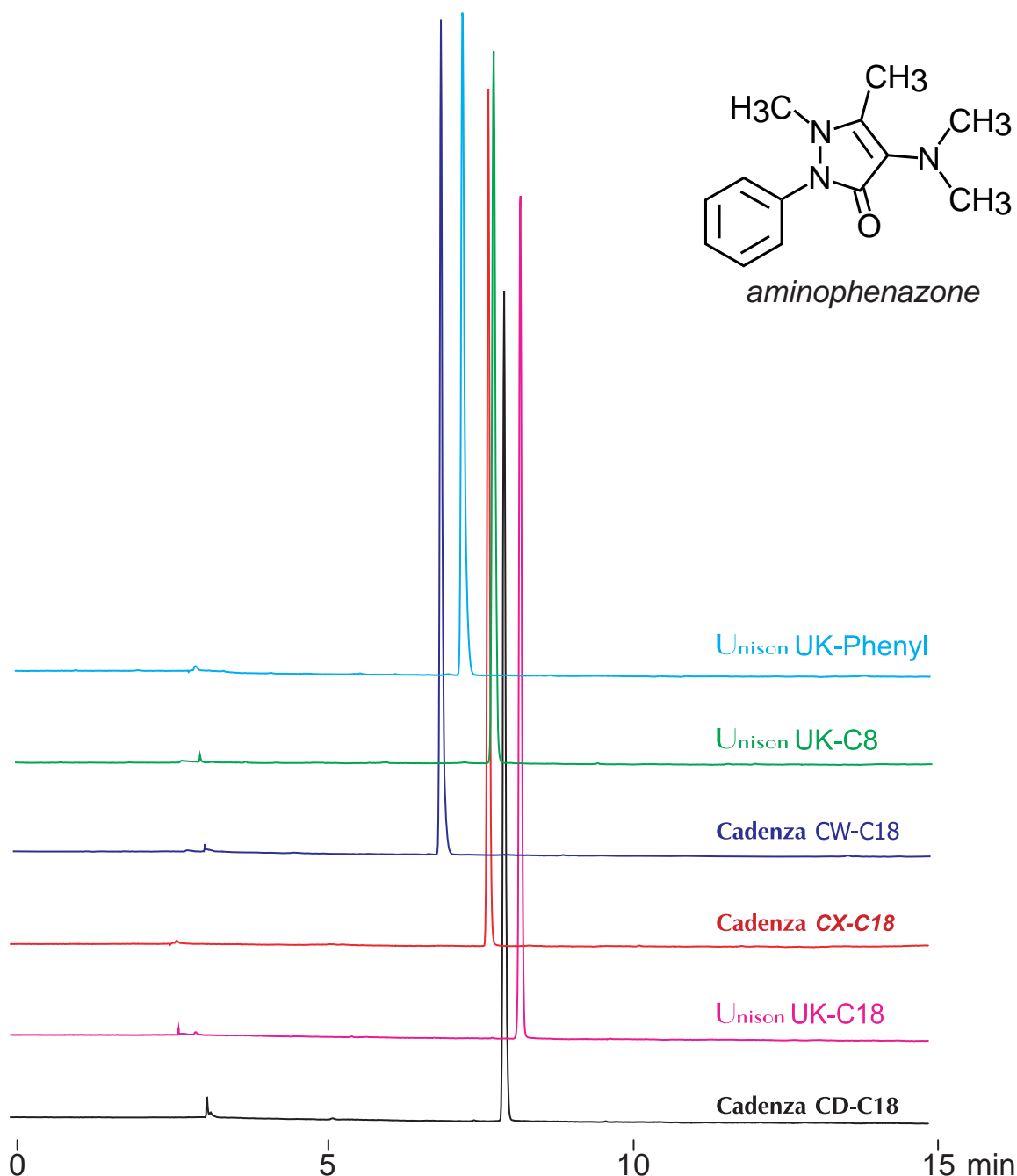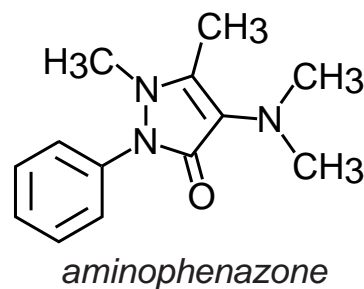

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-80%B (0-15min)

1.0mL/min (12MPa), 37°C, 260nm

10µL (0.1mg/mL, 100mM ammonium formate)