

Unison UK-Phenyl Cadenza CW-C18


Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Clarithromycin (Macrolide antibiotic)

クラリスロマイシン (マクロライド系抗生物質)

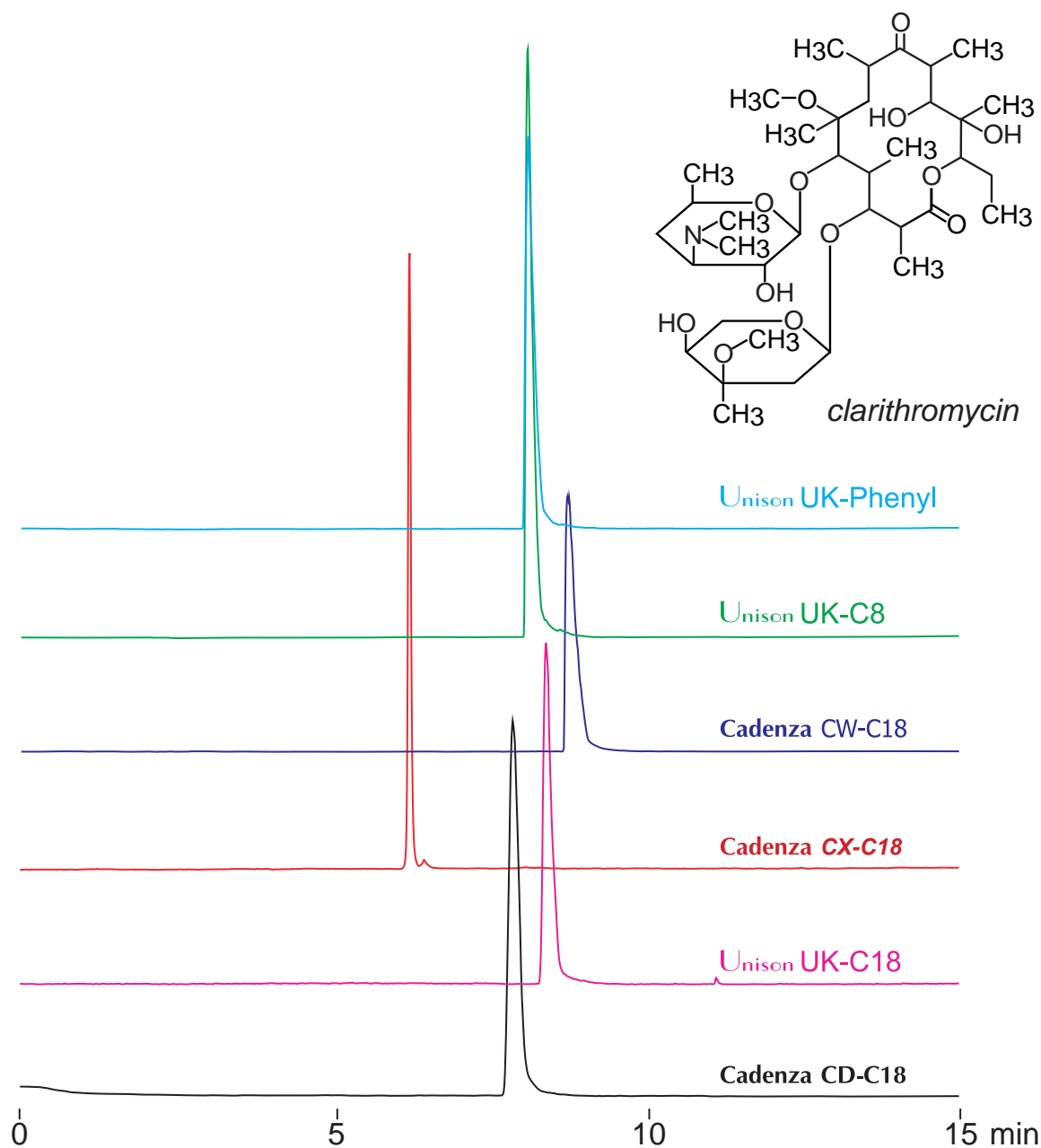

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-100%B (0-15min)

1.0mL/min (11MPa), 37°C

ELSD(spray chamber 30°C, drift tube 60°C)

10μL (1.0mg/mL, 100mM ammonium formate)