

Unison UK-Phenyl Cadenza CW-C18



Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Phenylbutazone (Anti-inflammatory)

フェニルブタゾン (ピラズロン誘導体, 消炎鎮痛薬)

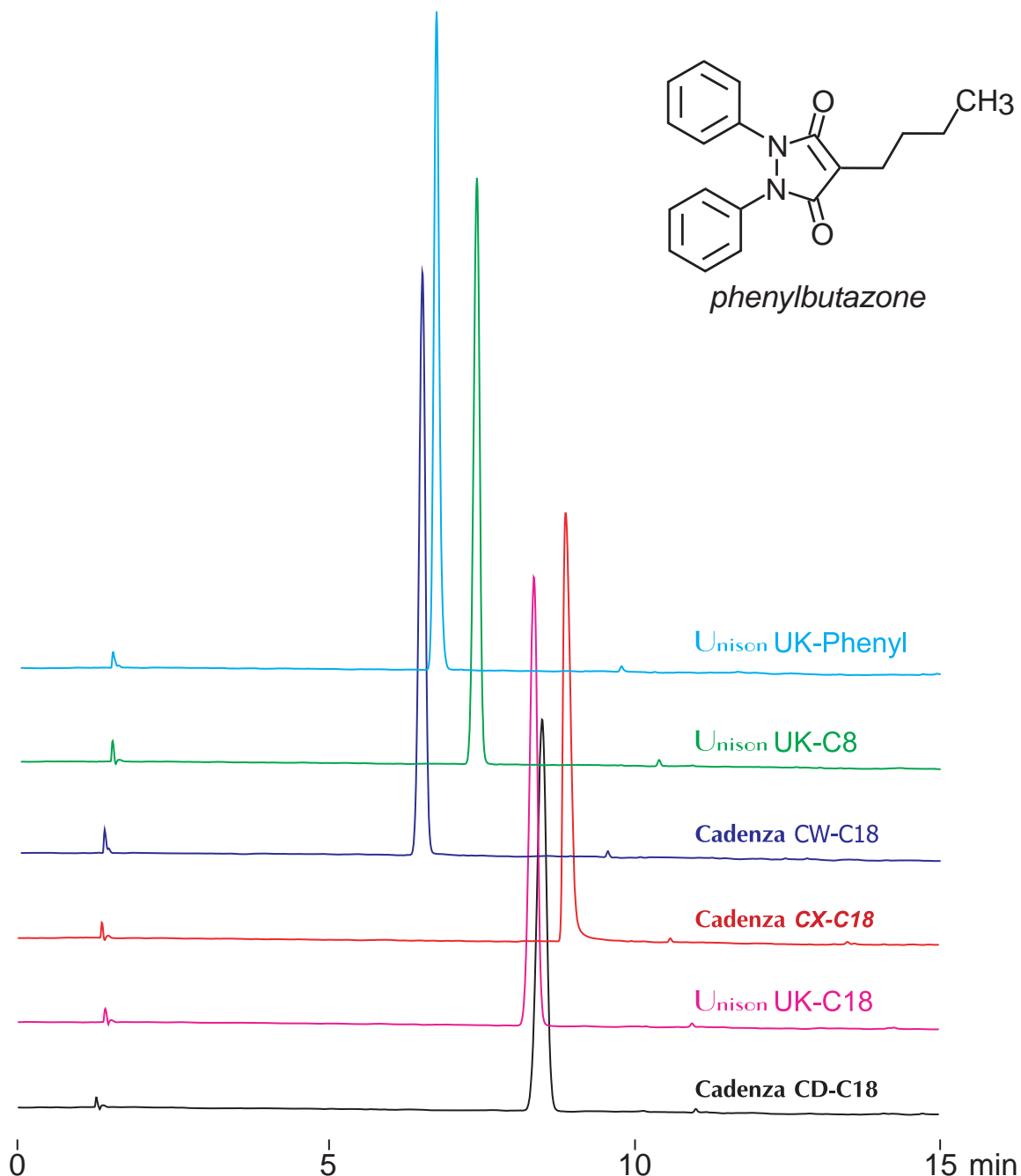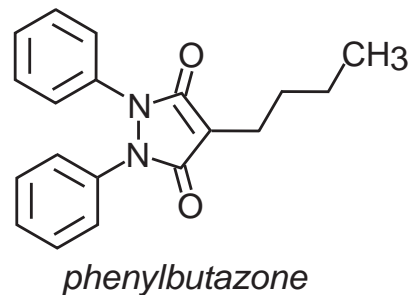

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

20-100%B (0-15min)

1.0mL/min (12MPa), 37°C, 250nm

10µL (0.1mg/mL, acetonitrile / 100mM ammonium formate = 50 / 50)