

Unison UK-Phenyl Cadenza CW-C18



Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Pitavastatin (Hypolipidemic agent)

ピタバスタチン (脂質異常症治療薬)

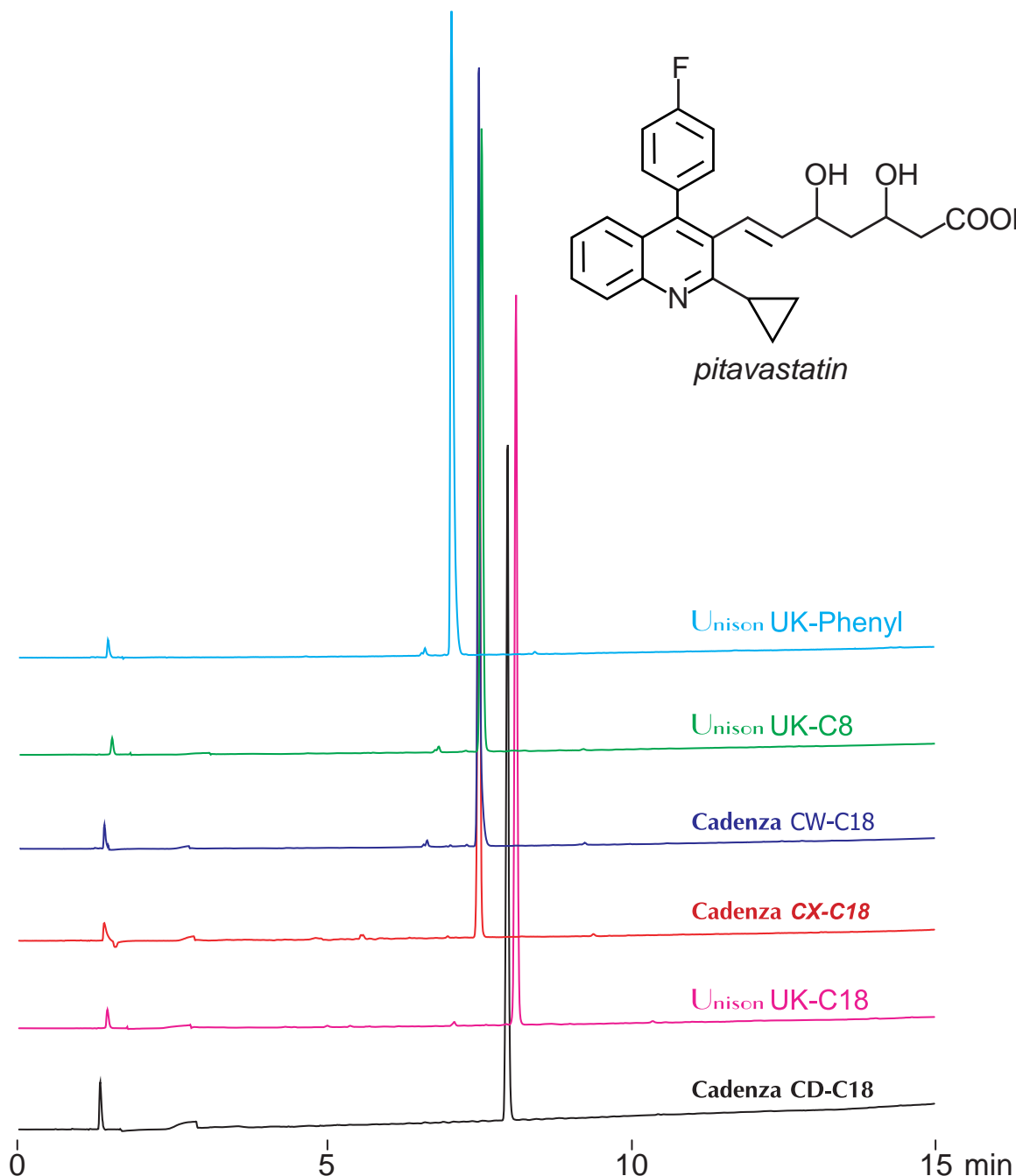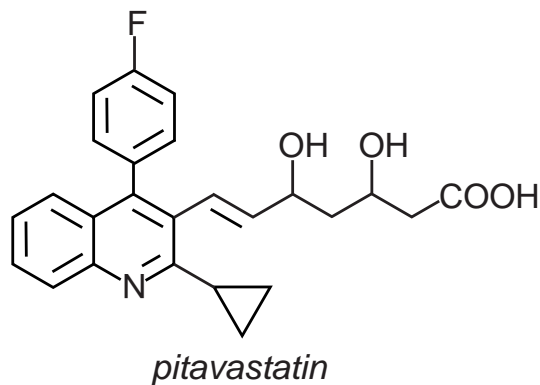

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-100%B (0-15min)

1.0mL/min (12MPa), 37°C, 250nm

10µL (0.1mg/mL, water / acetonitrile / formic acid = 50 / 50 / 1)