



Scherzo SW-C18
 Scherzo SM-C18
 Scherzo SS-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Chloronicotinic acid isomers

クロロニコチン酸異性体

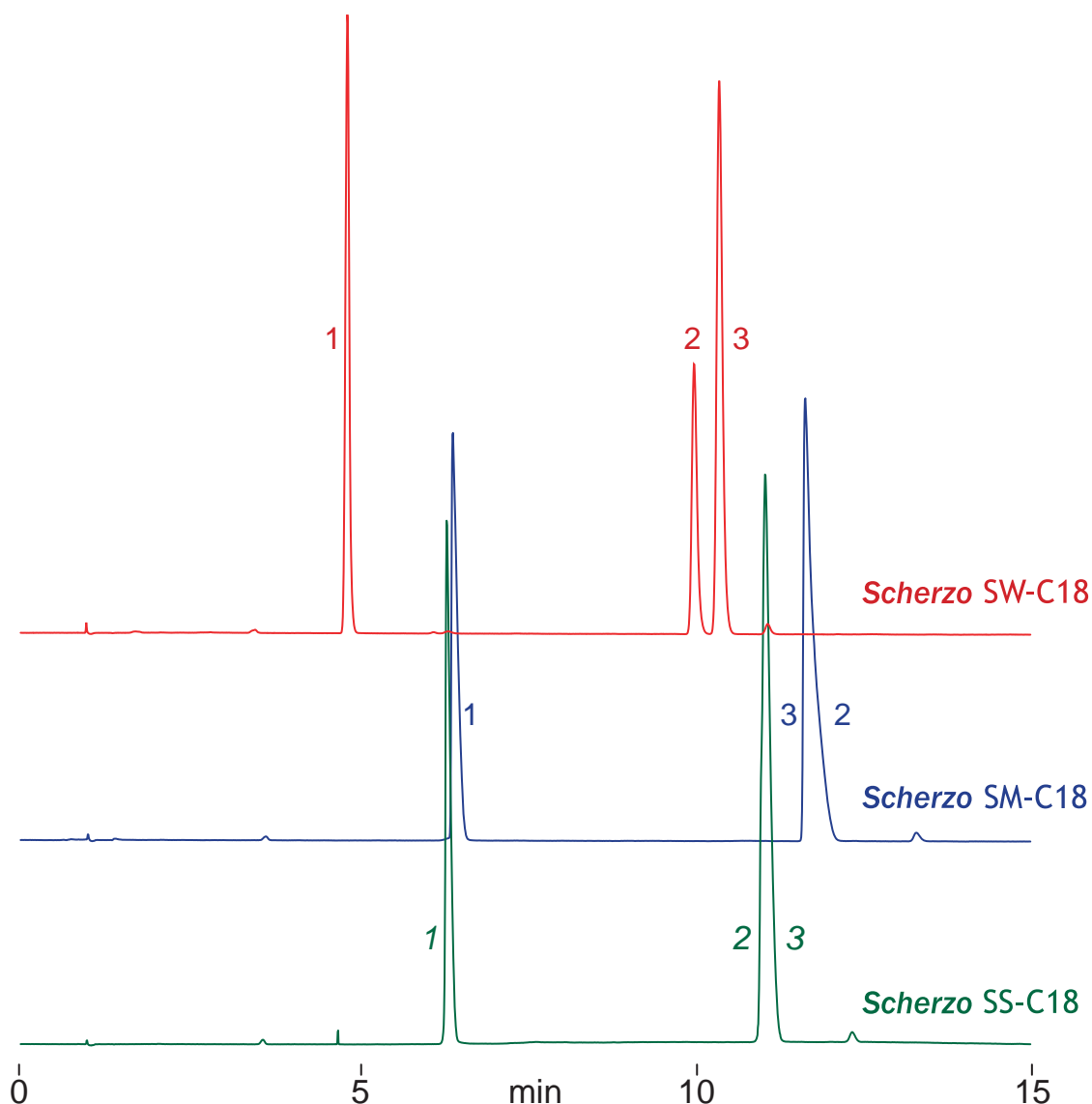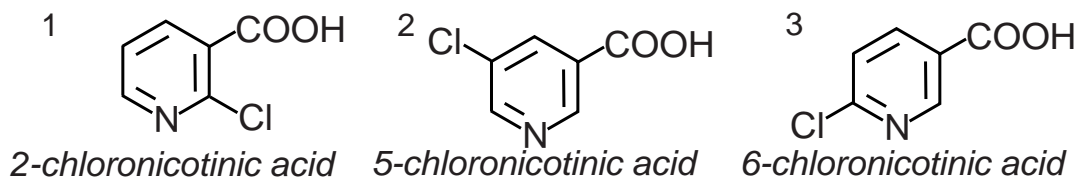

100 x 4.6 mm

A: water / formic acid = 100 / 0.3

B: 100mM ammonium formate / acetonitrile = 70 / 30

0-60%B (0-20min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.3mg/mL, water / methanol / formic acid = 50 / 50 / 1)