


Presto FF-C18

150 x 4.6 mm

HPLC COLUMN APPLICATION

Humic acid

フミン酸

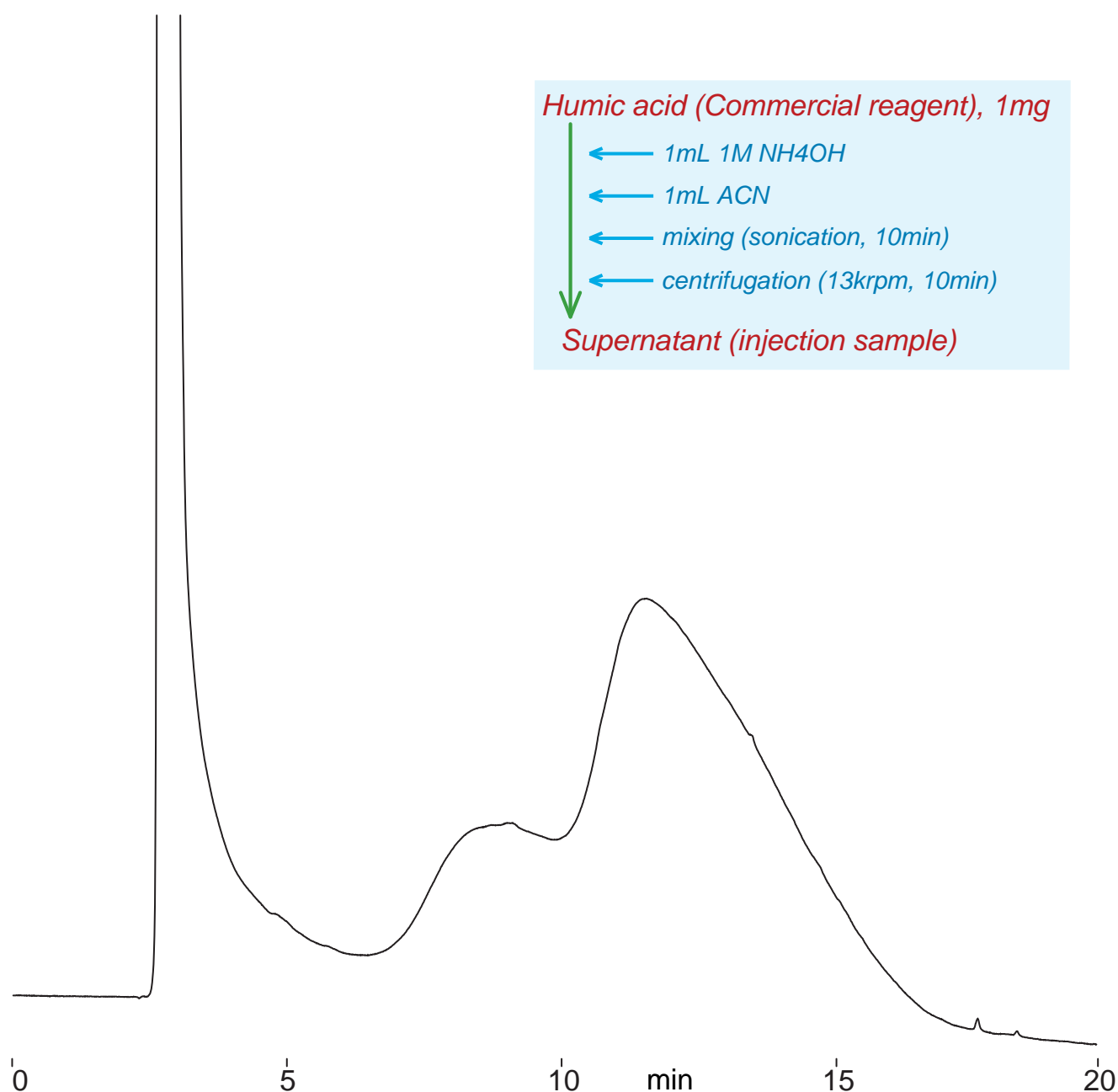

Presto FF-C18, 150 x 4.6 mm

A: 10mM AcONH₄B: ACN / 10mM AcONH₄ / 10mM NH₄OH = 90 / 9 / 1

20%B(0-1min), 20-100%B(1-15min), 100%B(15-20min)

0.4mL/min (21MPa), 37°C, 260nm

20μL (0.5mg/mL, 1M NH₄OH / ACN = 50 / 50)