



Scherzo SW-C18
 Scherzo SM-C18
 Scherzo SS-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Pinaverium (Antispasmodic)

ピナベリウム (鎮痙薬)

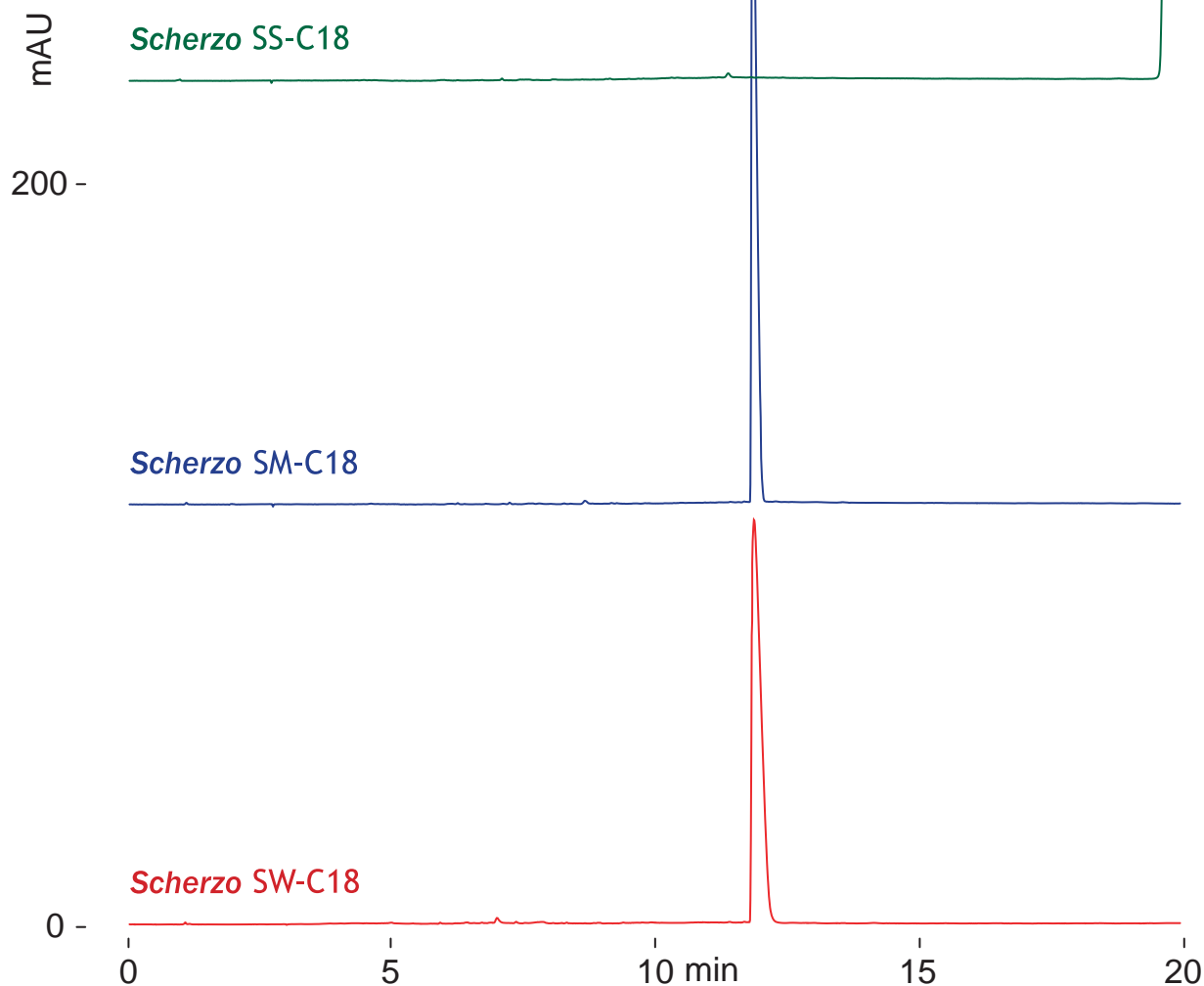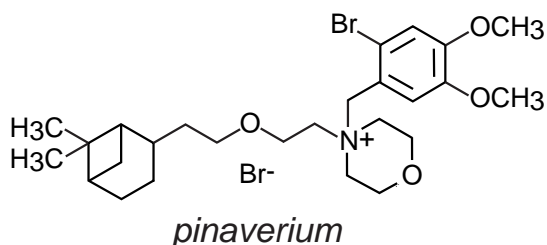

100 x 4.6 mm

A: 10mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 30 / 70

0-100%B (0-10min), 100%B (10-20min)

1.0mL/min (12MPa), 37°C, 250nm

10µL (1mg/mL, 100mM ammonium formate)