


Scherzo SM-C18

100 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Nisin, Natamycin (Food preservatives)

ナイシン, ナタマイシン (食品保存料) の LC-MS 分析

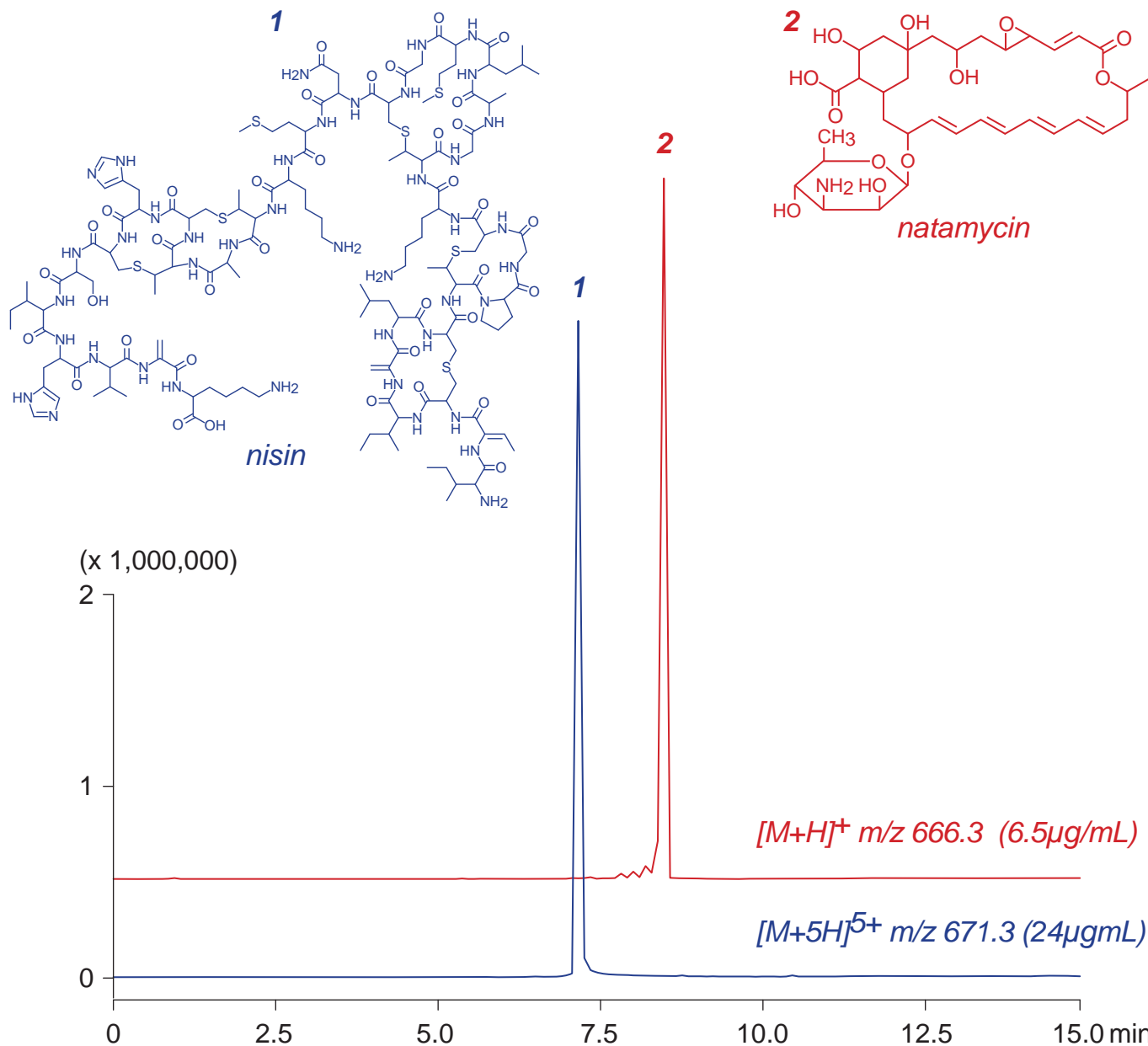

Scherzo SM-C18, 100 x 3 mm

A: water / formic acid = 100 / 0.5

B: acetonitrile / formic acid = 100 / 0.5

0-100%B (0-10min), 100%B (10-15min)

0.4mL/min (7MPa), 37°C, Single Quad. MS (ESI)

5µL (6.5-24µg/mL, 1% HCOOH)