

Cadenza **CX-C18** Unison **UK-Phenyl**

Unison **UK-C18** Unison **UK-C8**


Cadenza **CD-C18** Cadenza **CW-C18**

100 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Cyclosporin A (Antibiotics)

シクロスポリン A (抗生物質) の LC-MS 分析

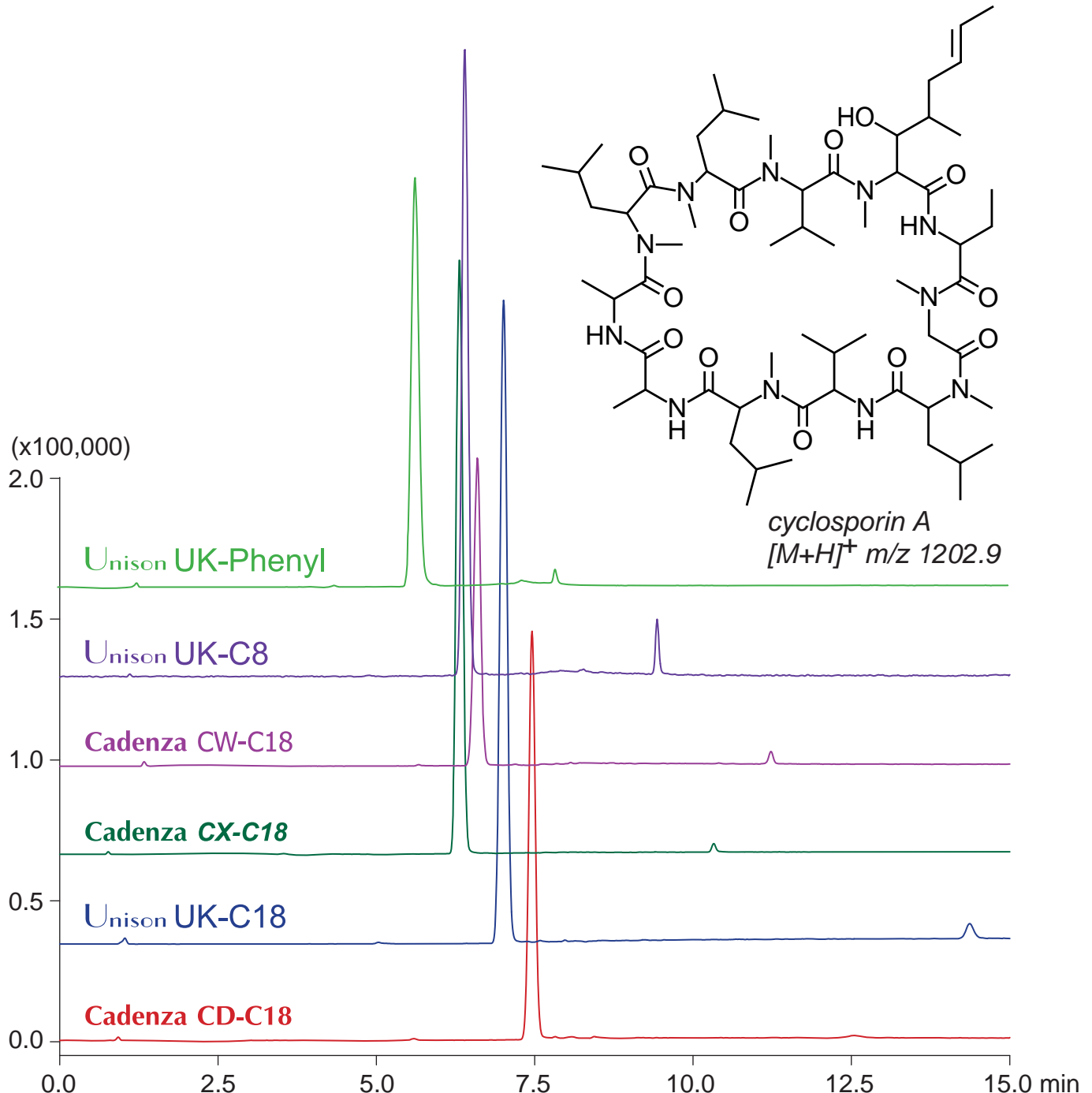

100 x 3 mm

A: water / formic acid = 100 / 0.1

B: methanol / formic acid = 100 / 0.1

60-100%B (0-6min), 100%B (6-15min)

0.4mL/min (6-9MPa), 50°C, Single Quad. MS (ESI)

5μL (2μg/mL, MeOH / water / HCOOH = 50 / 50 / 0.1)