

Unison UK-Phenyl Cadenza CW-C18


Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

N,N'-Diphenylethylenediamine (nitrogen-containing ligand)

N,N'-ジフェニルエチレンジアミン (含窒素配位子)

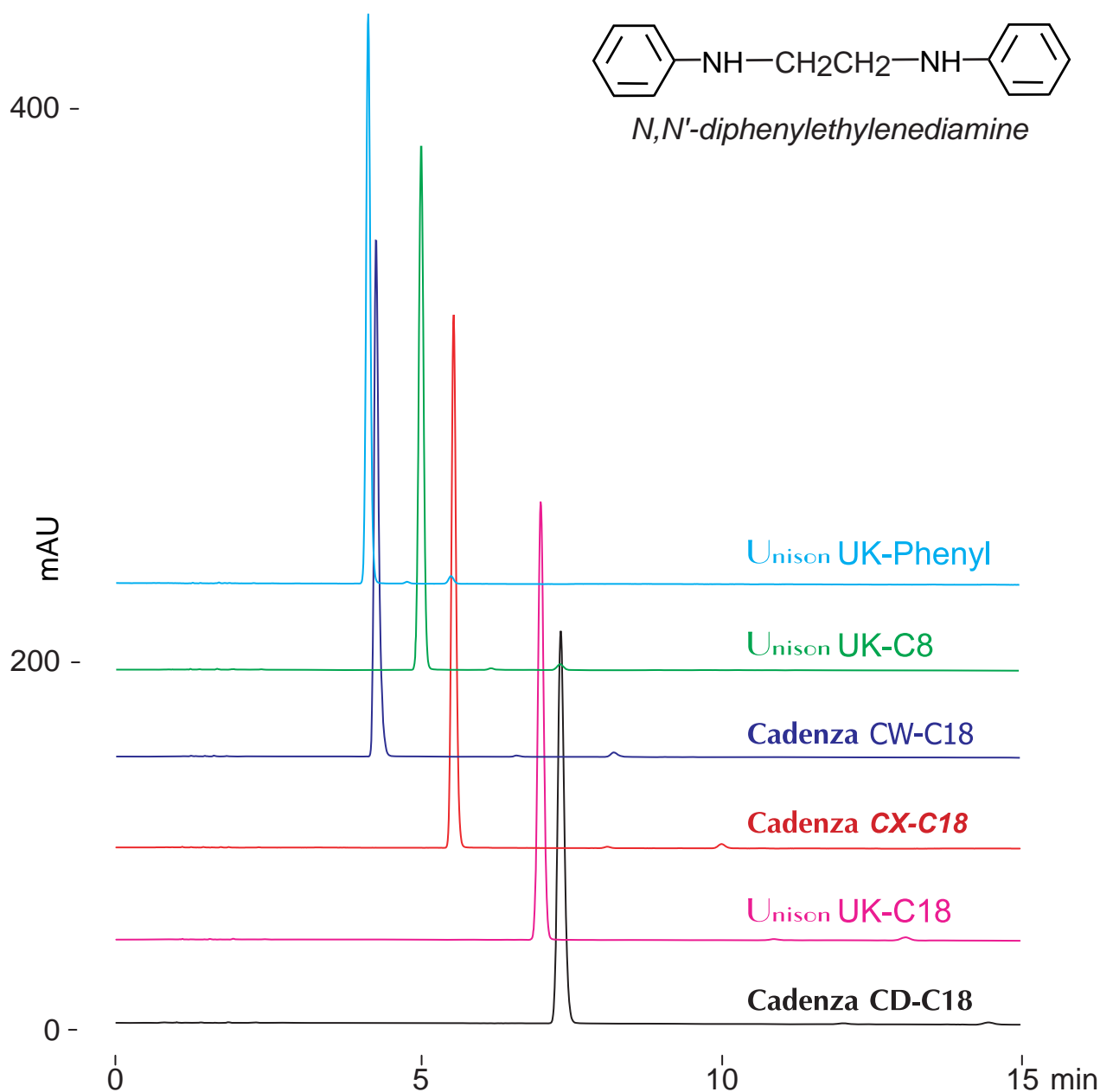

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 50 / 50

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-100%B (0-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.1mg/mL, 100mM ammonium formate / acetonitrile = 50 / 50)