

Unison UK-C8
 Cadenza CL-C18
 Cadenza CX-C18
 Cadenza CD-C18


Unison UK-Phenyl
 DACAPO DX-C18
 Cadenza CW-C18
 Unison UK-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Diosmin, Hesperidin (vascular protectant)

ジオスミン, ヘスペリジン (血管保護薬)

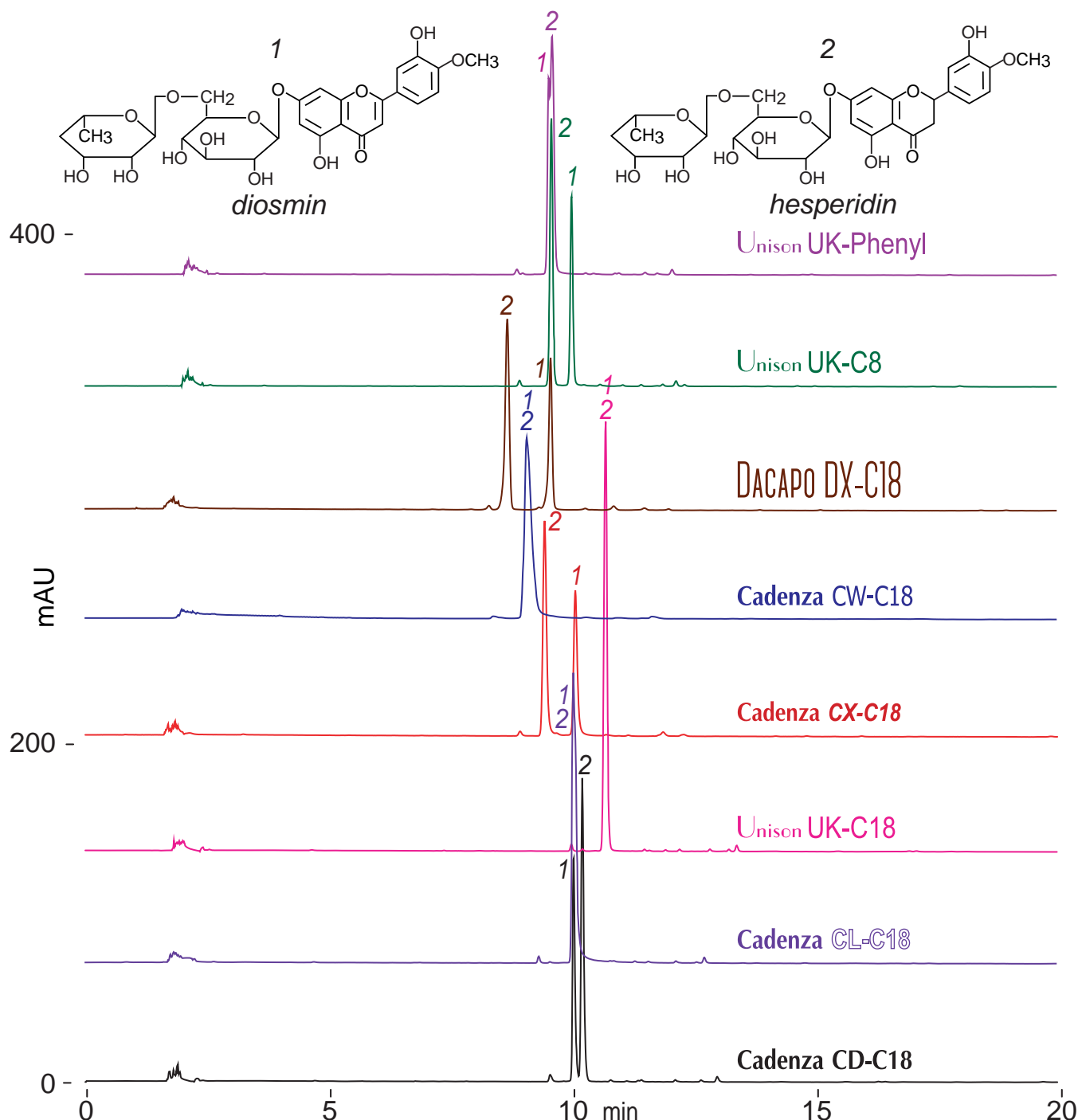

150 x 4.6 mm

A: water / acetic acid = 100 / 0.1

B: water / acetonitrile / acetic acid = 10 / 90 / 0.1

10-60%B (0-20min)

1.0mL/min (15-21MPa), 37°C, 280nm

10µL (0.125mg/mL, DMSO / MeOH / water / acetic acid = 50 / 25 / 25 / 1)