

Unison UK-C18

Cadenza CL-C18


Cadenza CD-C18

150 x 4.6 mm

HPLC COLUMN APPLICATION

Androstenedione biosynthetic pathway

アンドロステンジオン生合成経路

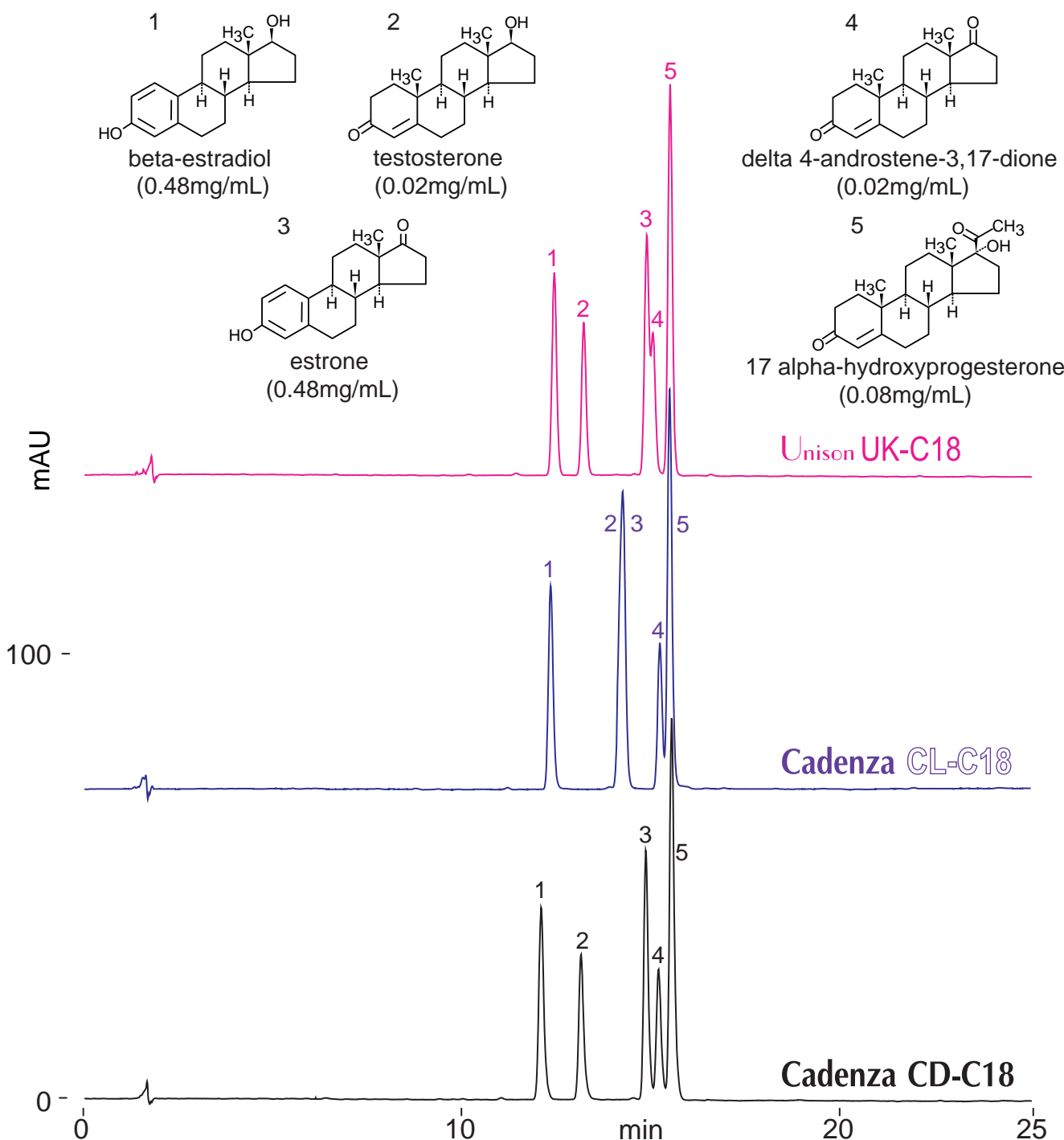

150 x 4.6 mm

A: water / acetic acid = 100 / 0.5

B: acetonitrile / acetic acid = 100 / 0.5

35-75%B (0-25min)

1.0mL/min (15MPa), 37°C, 260nm

10µL (0.02-0.48mg/mL, ACN / formic acid = 100 / 1)