


Scherzo SM-C18

150 x 4.6 mm

HPLC COLUMN APPLICATION

Quinolinol isomers

キノリノール異性体

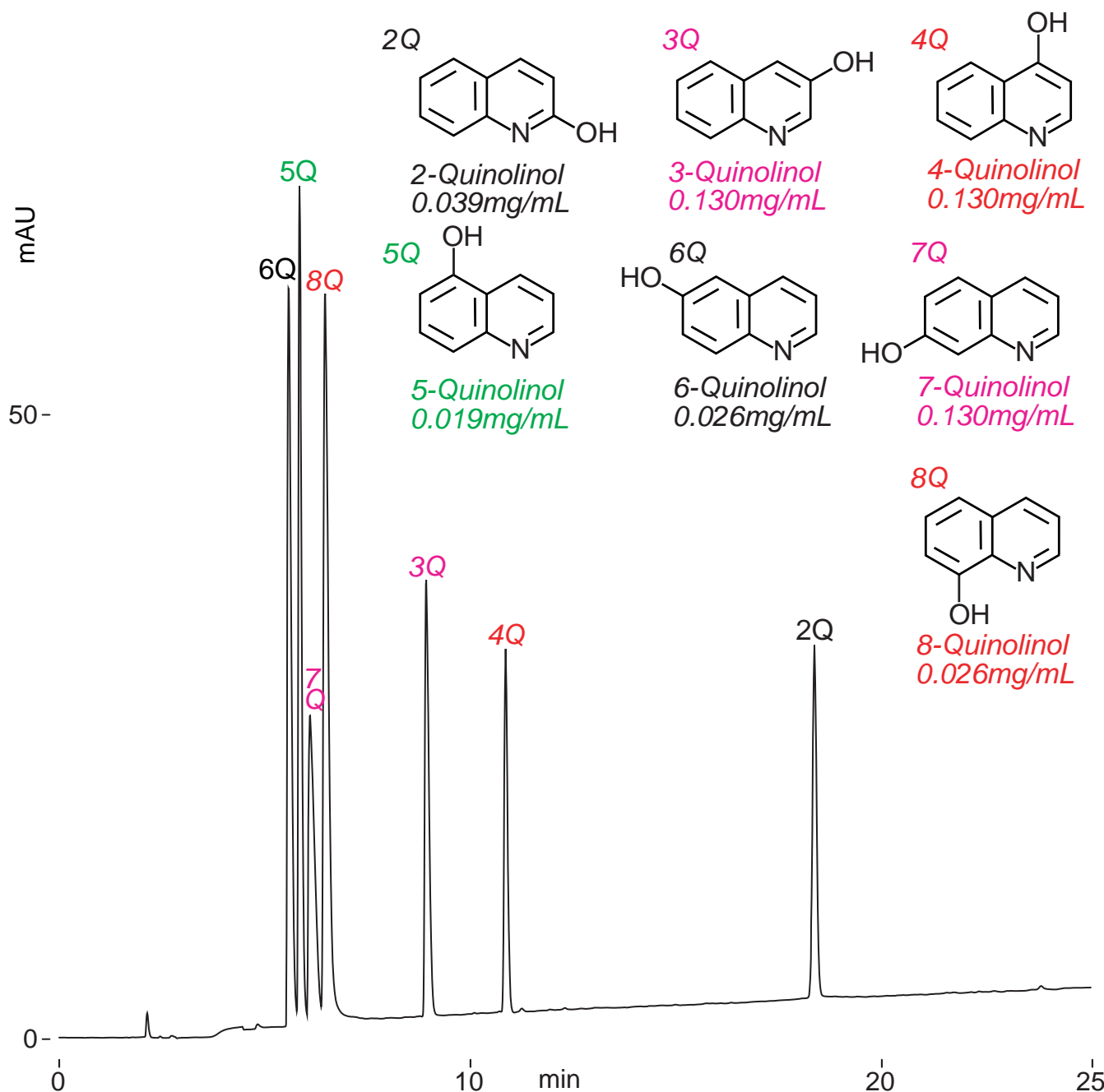

Scherzo SM-C18, 150 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 70 / 30 / 0.5

0-100%B (0-25min)

1.0mL/min (11MPa), 37°C, 260nm

10µL(0.026-0.130mg/mL, water / ACN / HCOOH = 74 / 3 / 1)