

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18

Cadenza CW-C18


Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Entacapone (Antiparkinsonian)

エンタカポン (パーキンソン病治療薬)

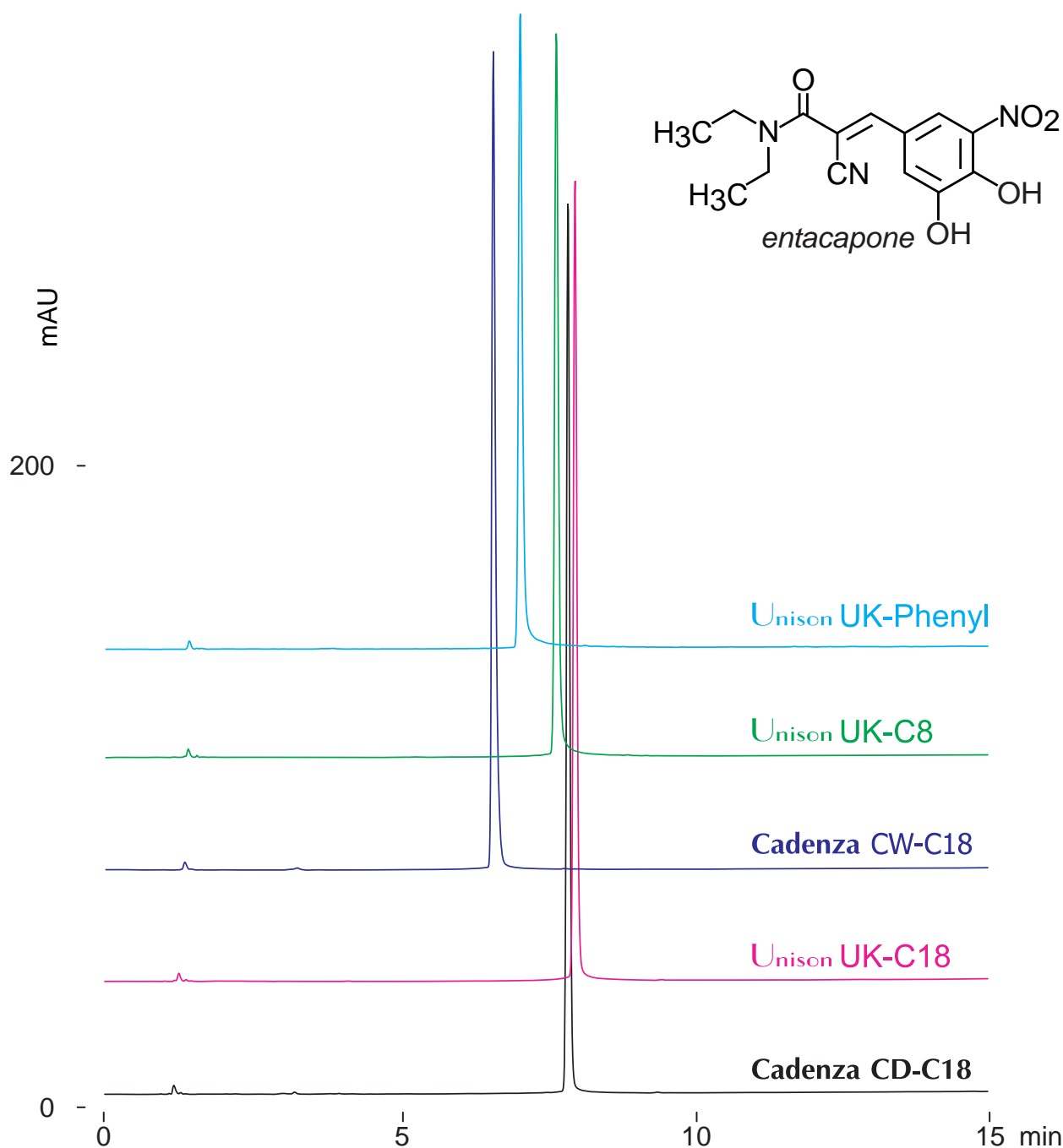

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

20-100%B (0-15min)

1.0mL/min (11MPa), 37°C, 260nm

10µL (0.25mg/mL, acetonitrile / water / formic acid = 50 / 50 / 1)