

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18



Cadenza CW-C18

Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Loxoprofen (Nonsteroidal anti-inflammatory drug)

ロキソプロフェン (非ステロイド性抗炎症薬)

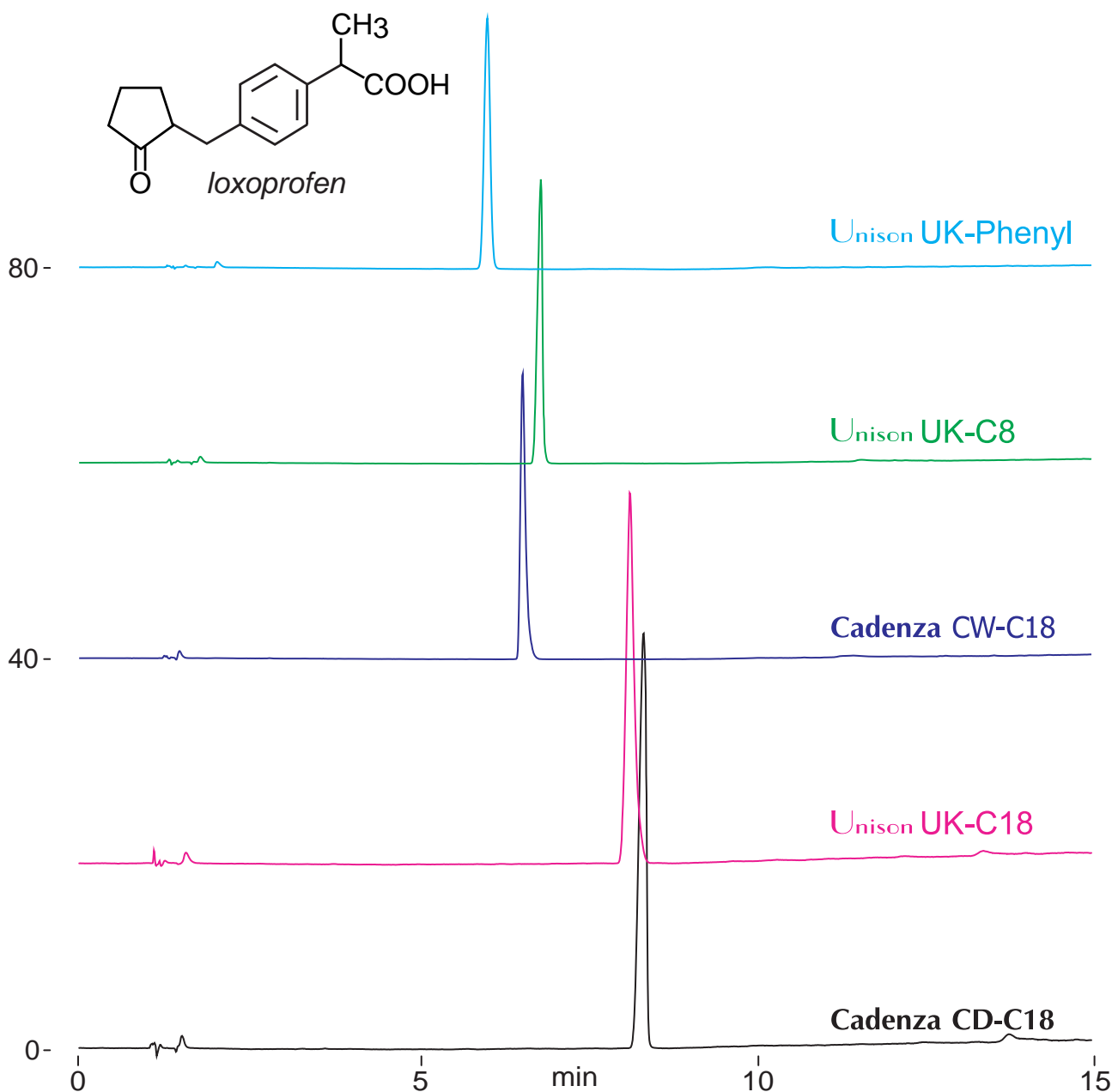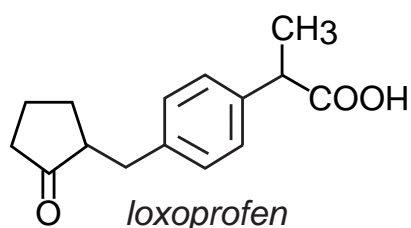

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / methanol = 10 / 90

25-100%B (0-15min,)

1.0mL/min (16MPa), 37°C, 260nm

10µL (1.0mg/mL, 100mM ammonium formate)