

Unison UK-Phenyl Cadenza CW-C18



Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Suplatast (Antiallergic)

スプラタスト (抗アレルギー薬)

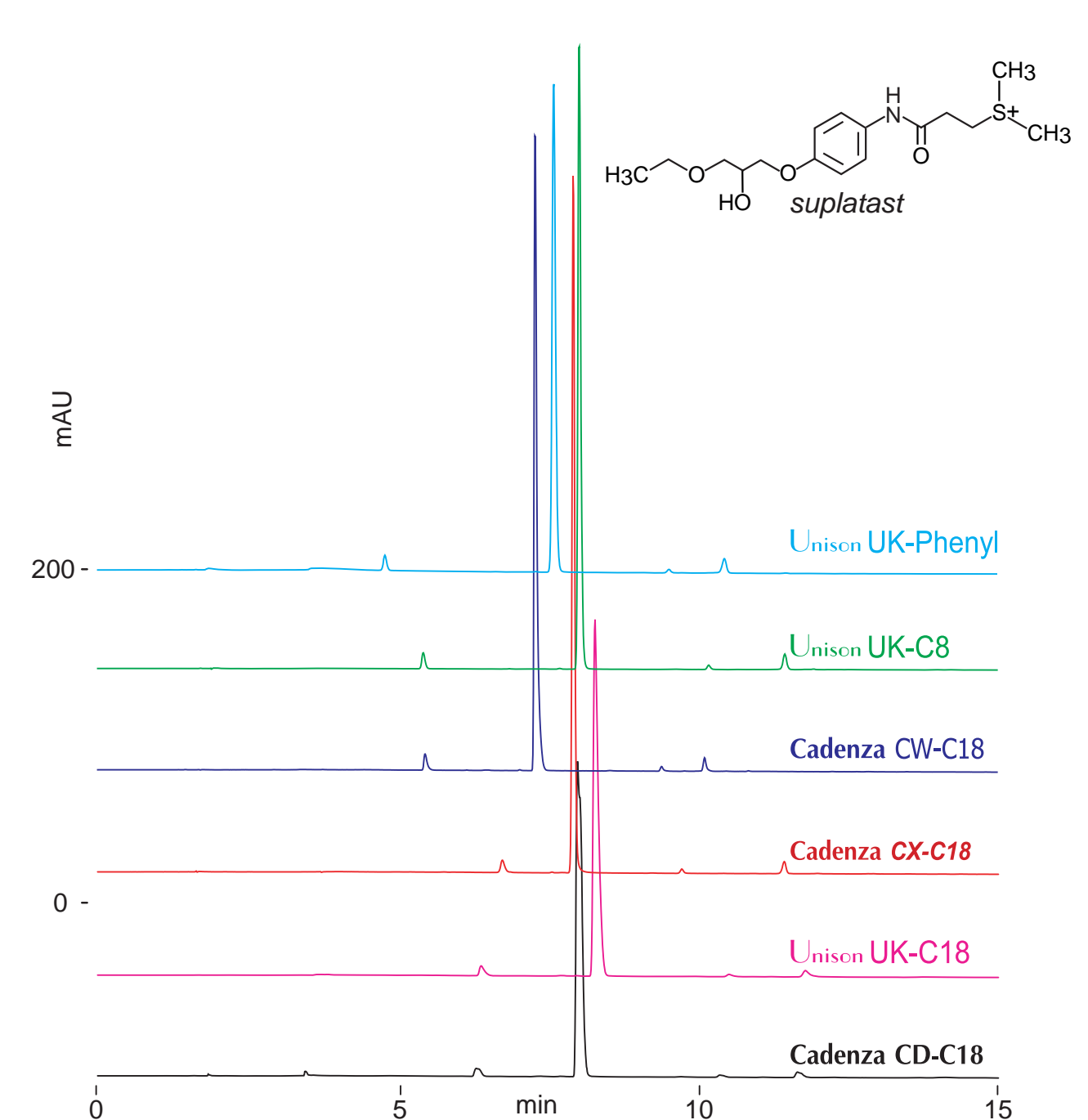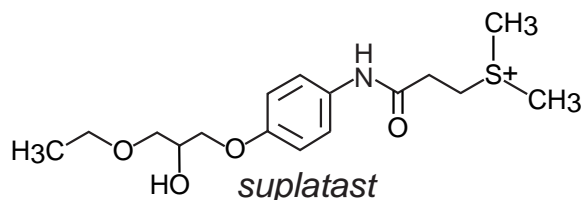

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-50%B (0-15min)

1.0mL/min (14MPa), 37°C, 260nm

10µL (1.0mg/mL, 100mM ammonium formate)