

Unison UK-Phenyl Cadenza CW-C18


Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

### Quinoline, isoquinoline

キノリン, イソキノリン

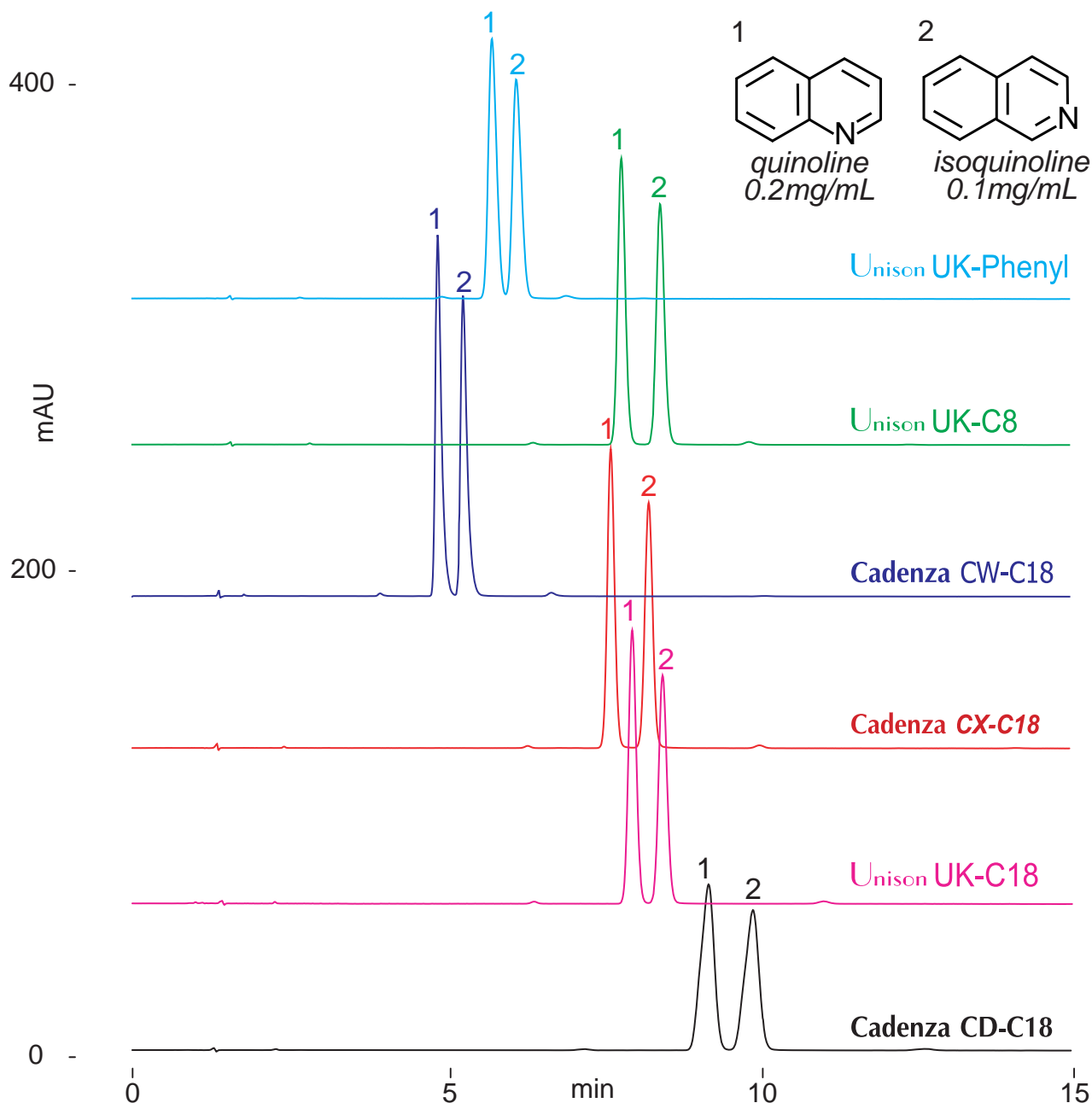

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

25-35%B (0-15min)

1.0mL/min (15MPa), 37°C, 260nm

10µL (0.1-0.2mg/mL, 100mM ammonium formate)