




Scherzo SW-C18  
 Scherzo SM-C18  
 Scherzo SS-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Adrenaline, lidocain (Local anesthetics)

アドレナリン, リドカイン (局所麻酔薬)

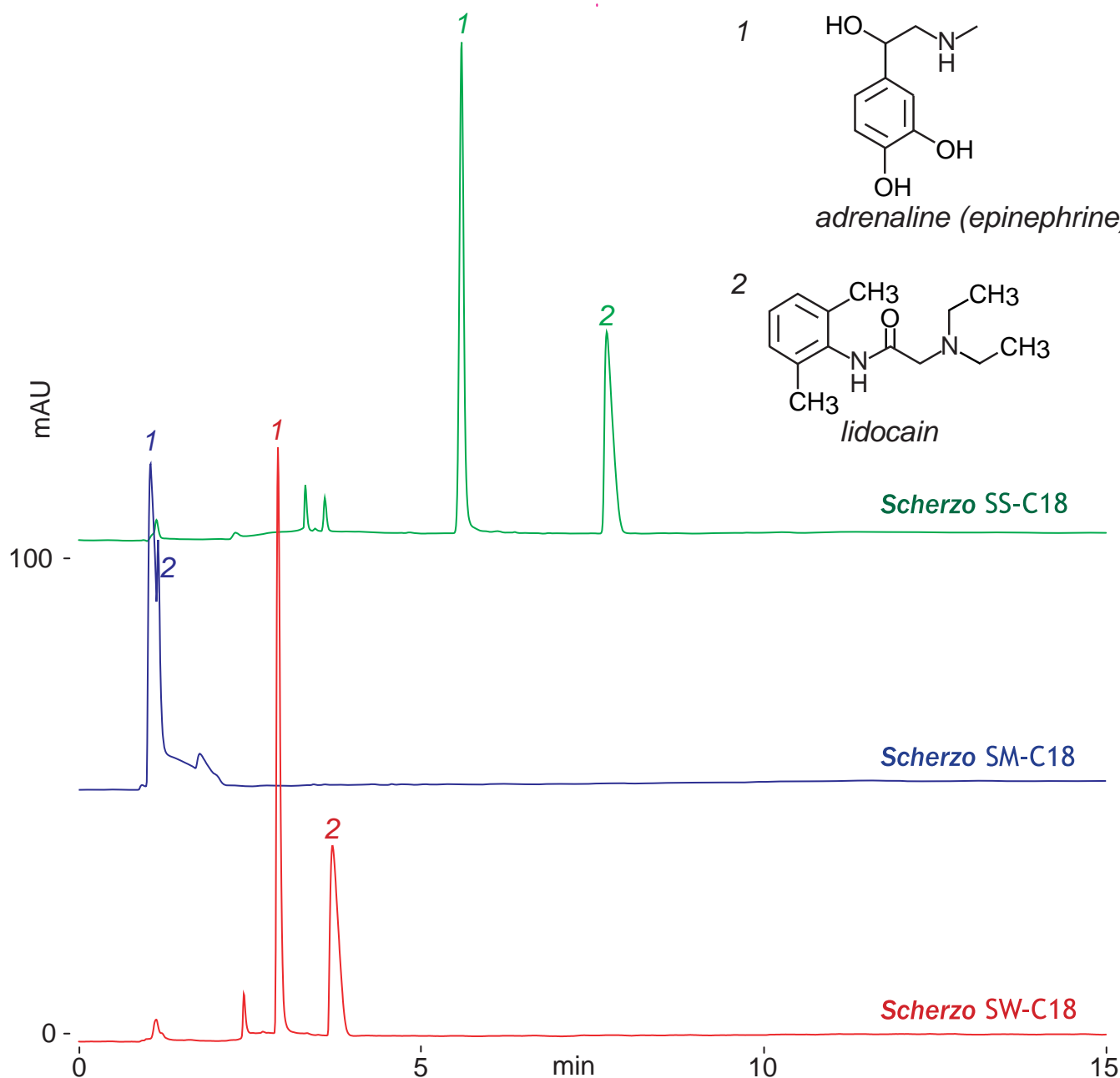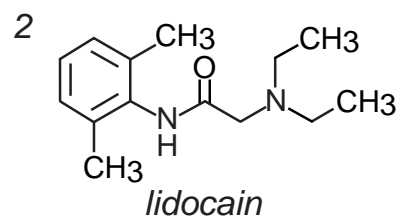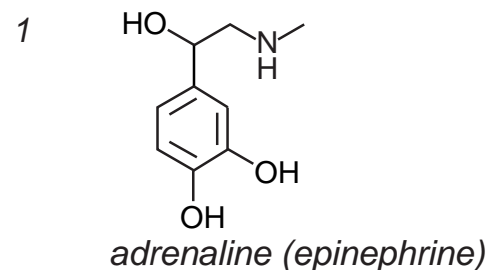

100 x 4.6 mm

A: water / acetonitrile / acetic acid = 50 / 50 / 0.1

B: 100mM ammonium acetate / acetonitrile = 30 / 70

0-100%B (0-10min), 100%B(10-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.5mg/mL, water / acetic acid = 100 / 1)