


Scherzo SS-C18

100 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Adrenaline, lidocaine (local anesthetics)

アドレナリン, リドカインの (局所麻酔薬) LC-MS 分析

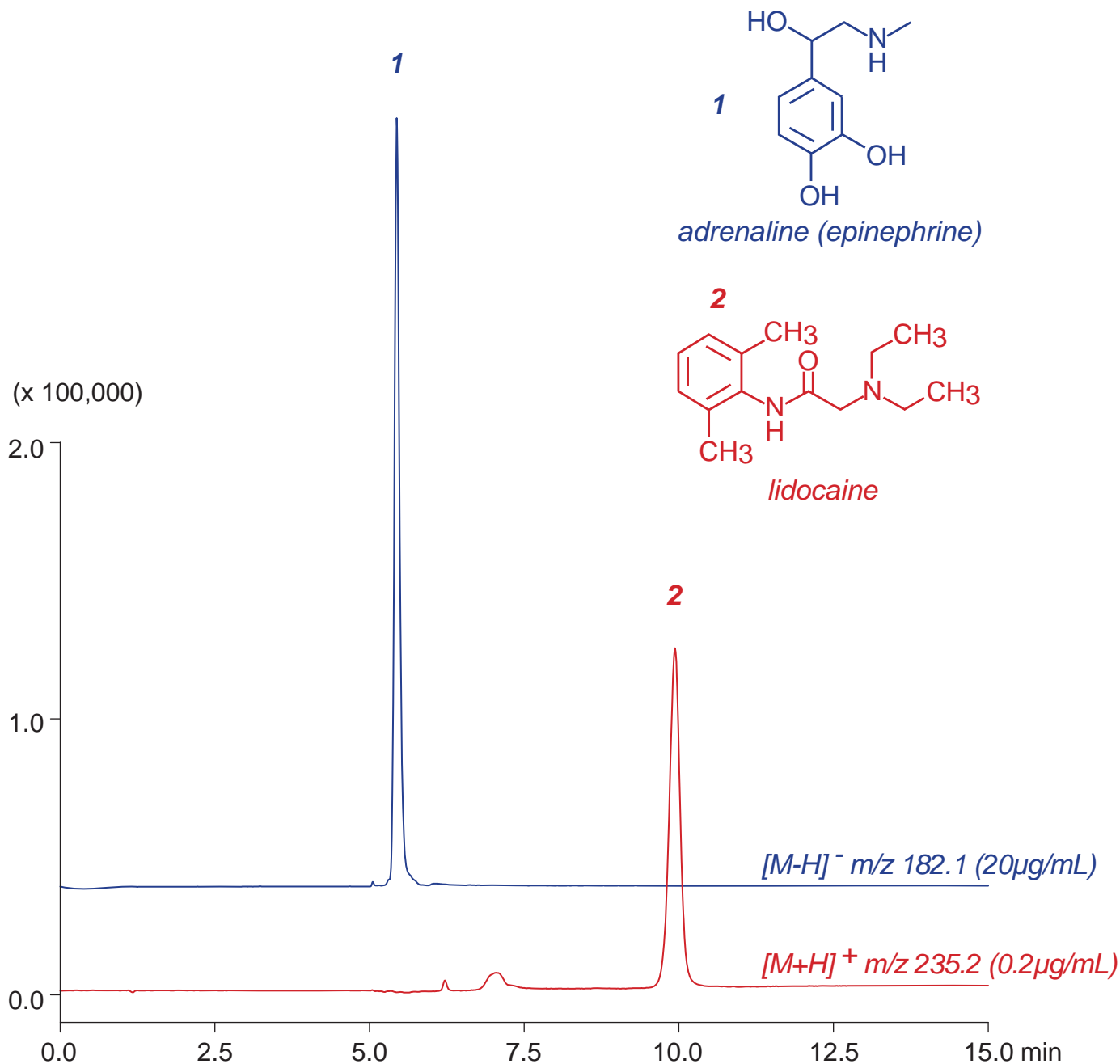

Scherzo SS-C18, 100 x 3 mm

A: methanol / water / formic acid = 50 / 50 / 0.1

B: methanol / 100mM ammonium formate = 50 / 50

0-100%B (0-10min), 100%B (10-15min)

0.4mL/min (11MPa), 37°C, Single Quad. MS (ESI)

5µL (0.2-20µg/mL, 0.5% HCOOH)