


Scherzo SM-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Adrenaline, lidocain (Local anesthetics)

アドレナリン, リドカイン (局所麻酔薬)

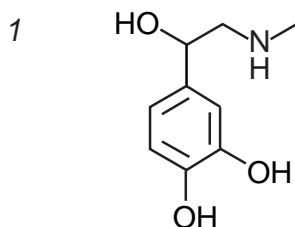


adrenaline (epinephrine)

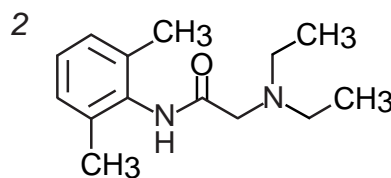

lidocain

100 -

mAU

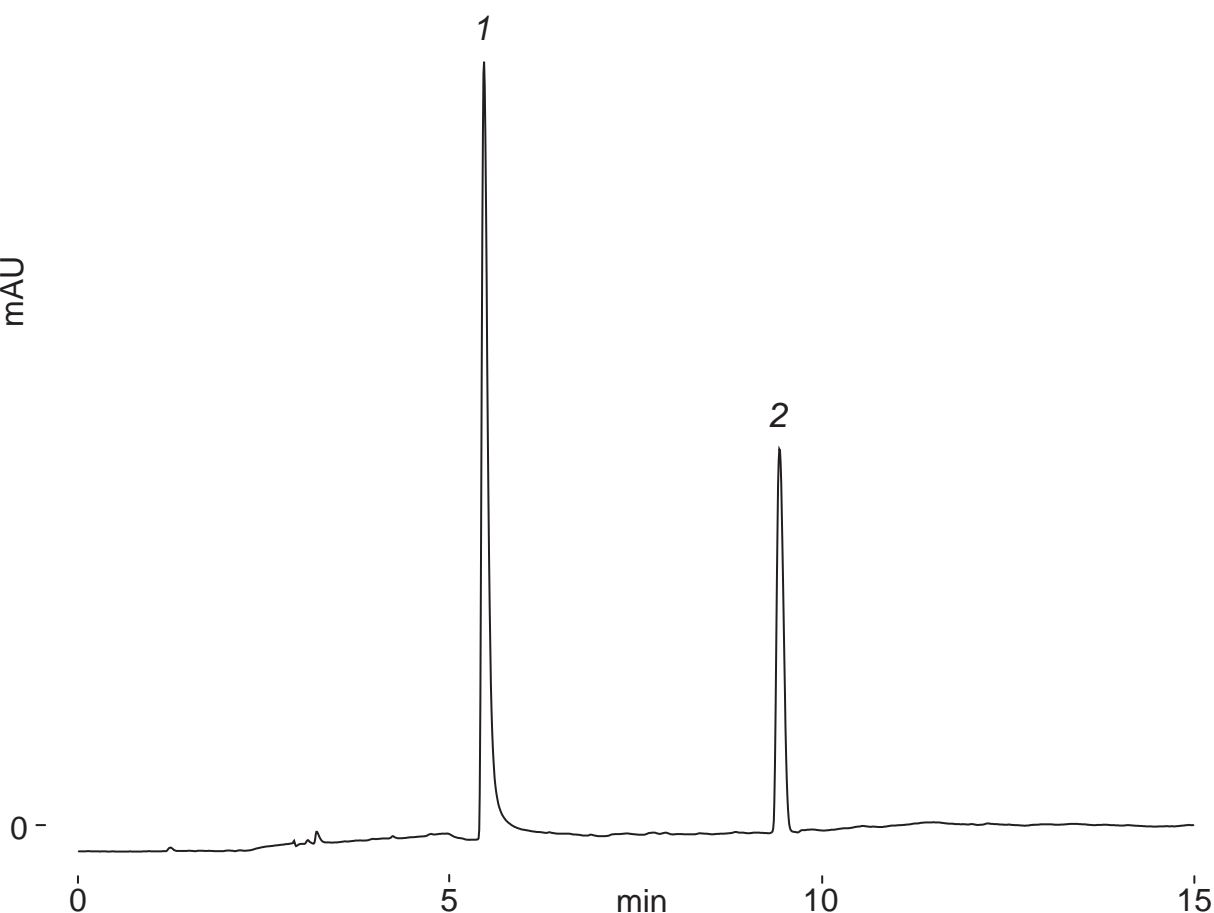

100 x 4.6 mm

A: 5mM ammonium acetate

B: 100mM ammonium acetate / acetonitrile = 30 / 70

0-100%B (0-10min), 100%B (10-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.5mg/mL, water / acetic acid = 100 / 0.1)