


Cadenza CX-C18

100 x 3 mm


HPLC COLUMN APPLICATION

LC-MS: Sibutramine (Anoretics)

シブトラミン (食欲抑制作用薬) の LC-MS 分析

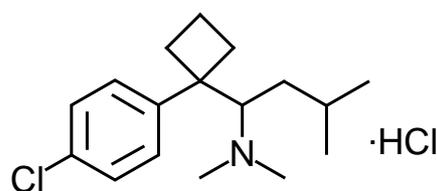

sibutramine hydrochloride

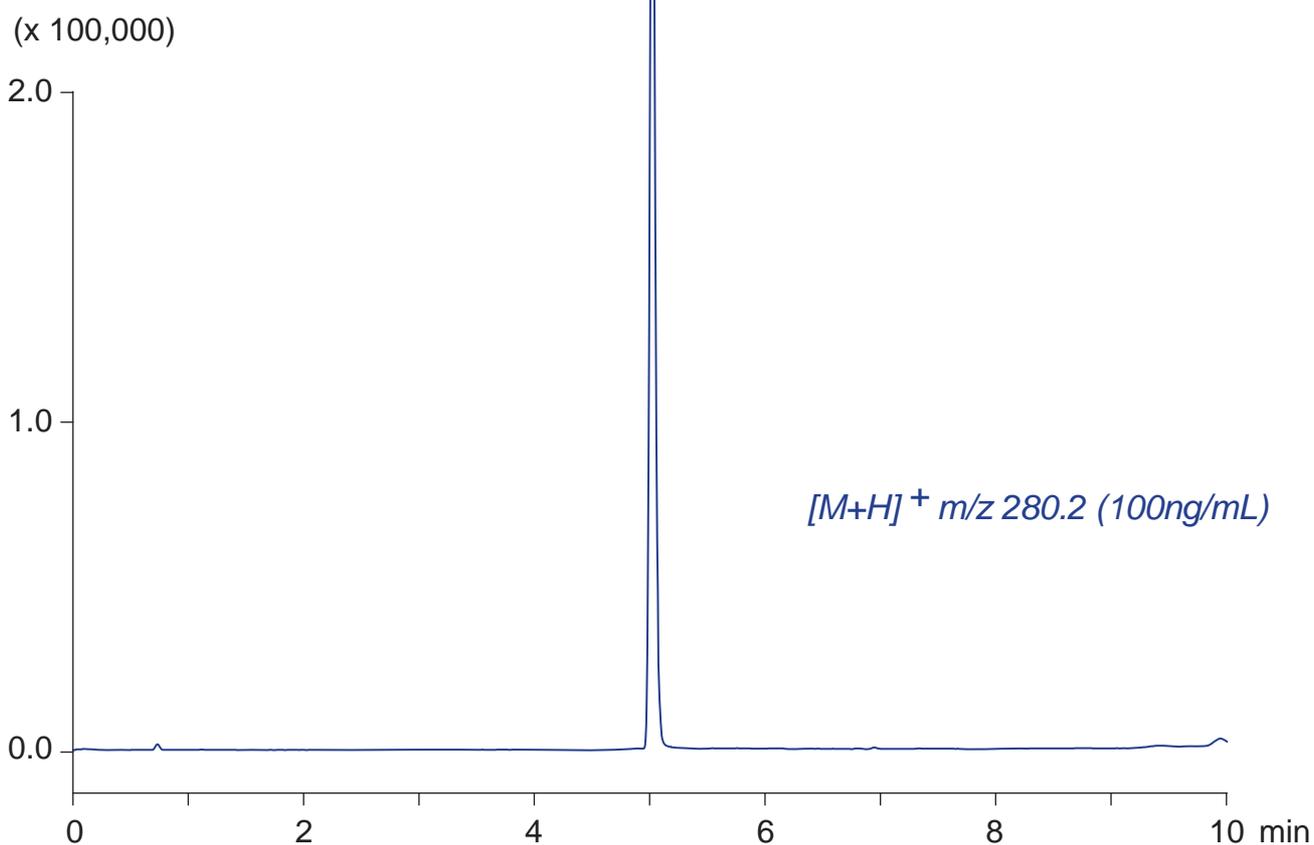

Cadenza CX-C18, 100 x 3 mm

A: acetonitrile / water / formic acid = 10 / 90 / 0.1

B: acetonitrile / water / formic acid = 90 / 10 / 0.1

0-60%B (2-10min)

0.4mL/min (7MPa), 37°C, Single Quad. MS (ESI)

5µL (100ng/mL, 0.5% HCOOH)