


DACAPO DX-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Sibutramine (Anoretics)

シブトラミン (食欲抑制作用薬)

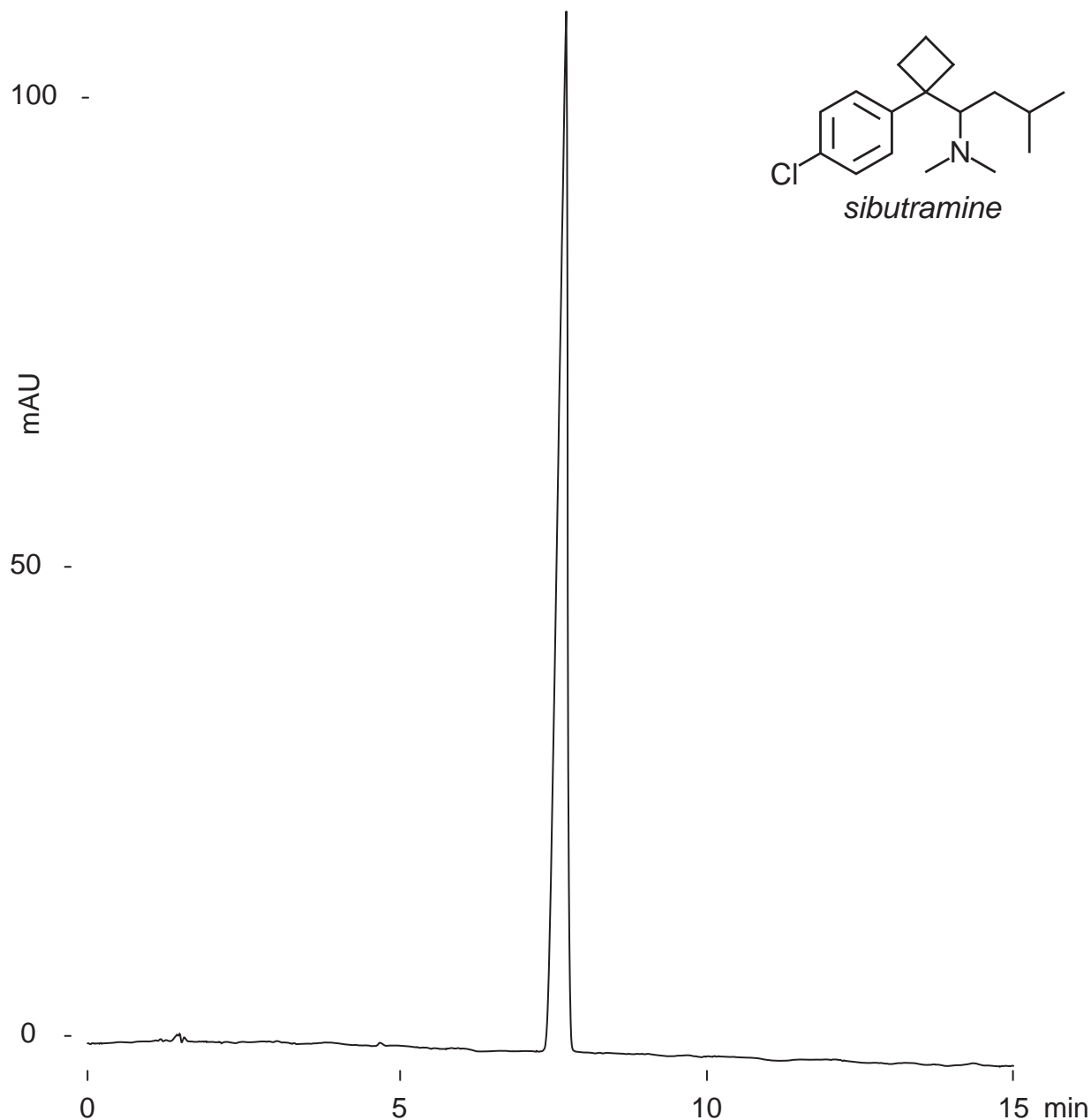

Dacapo DX-C18, 100 x 4.6 mm
 A: ACN / 10mM NH₄OH = 10 / 90
 B: ACN / 10mM NH₄OH = 90 / 10
 65-100%B (0-15min)
 1mL/min (12MPa), 37°C, 260nm
 10μL (1mg/mL, ACN / 10mM NH₄OH = 50 / 50)