


Scherzo SW-C18  
 Scherzo SM-C18  
 Scherzo SS-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Styrene urinary metabolites

尿中スチレン代謝物

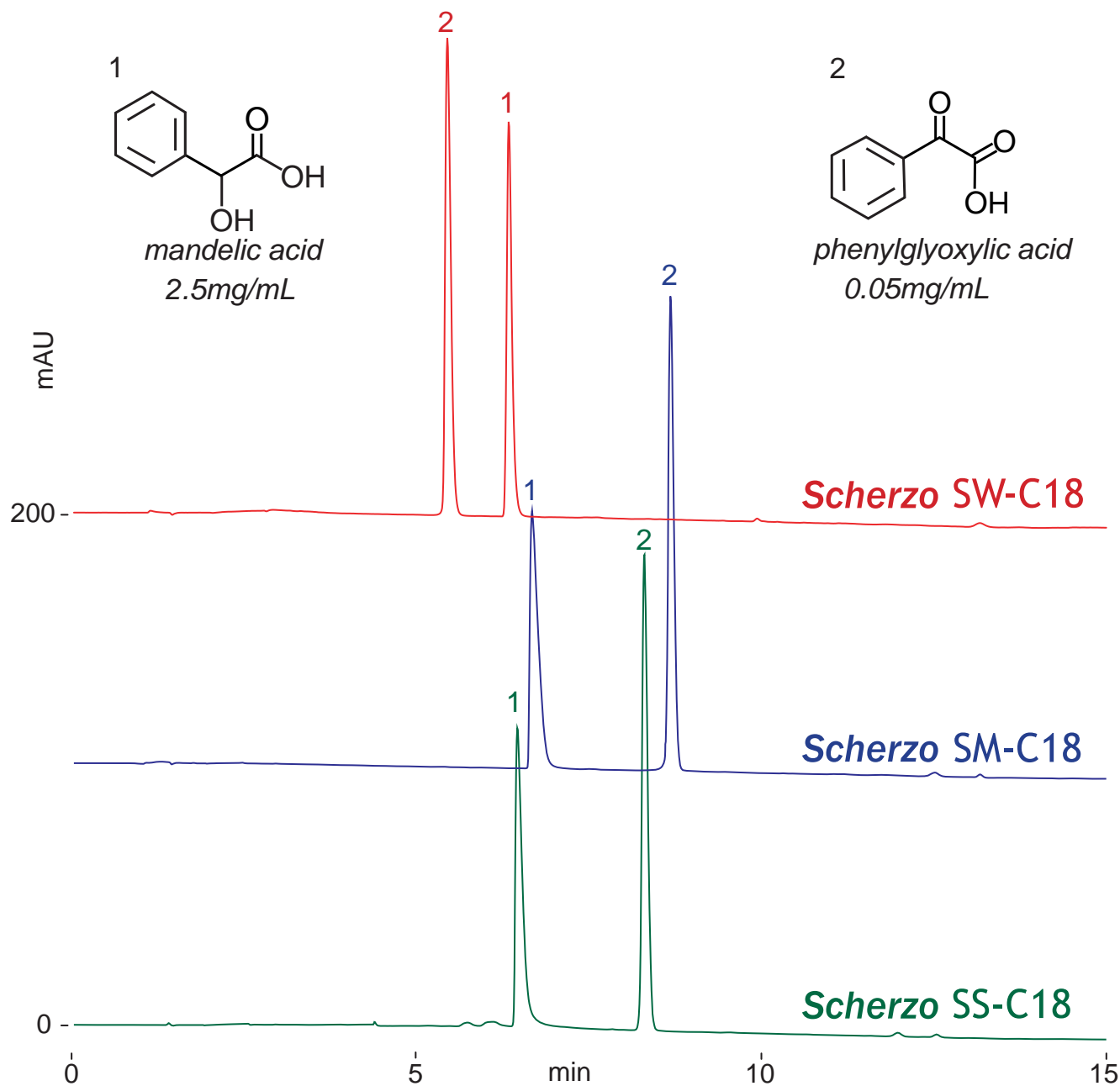

100 x 4.6 mm

A: water / HCOOH = 100 / 0.3

B: 100mM ammonium formate / acetonitrile = 70 / 30

0-100%B (0-12min), 100%B(12-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.05-2.5mg/mL, water / HCOOH = 100 / 1)