



Intrada Organic Acid

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS: Dobesilate (Vasoprotectives)

ドベシル酸 (血管保護薬) の LC-MS 分析

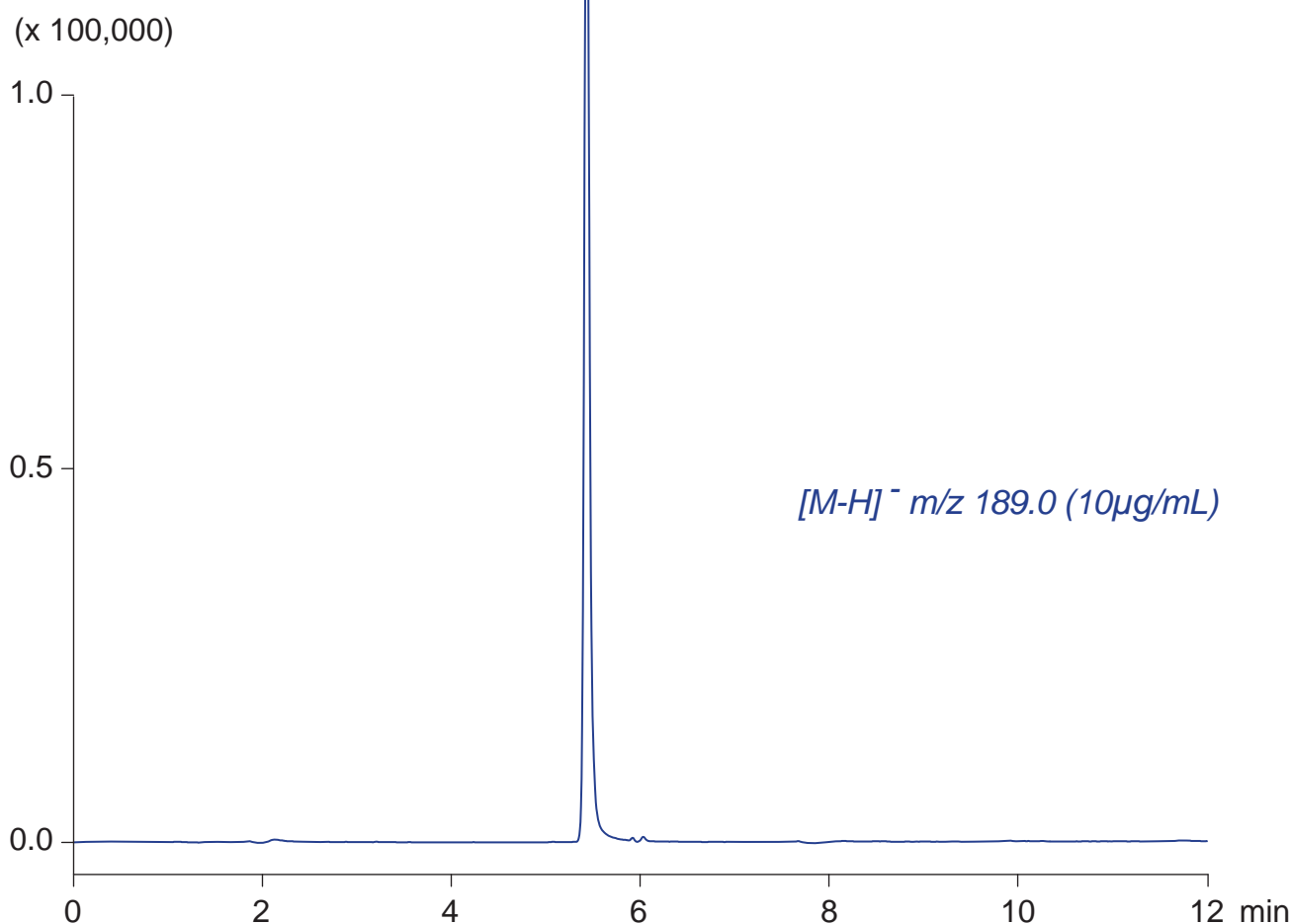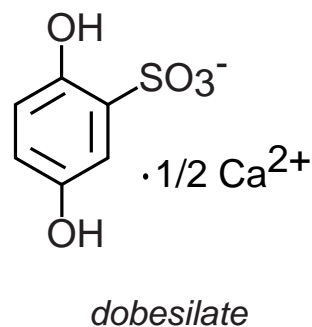

Intrada Organic Acid, 150 x 2 mm

A: acetonitrile / water / formic acid = 10 / 90 / 0.1

B: acetonitrile / 100mM ammonium formate = 10 / 90

0-100%B (0-6min), 100%B (6-12min)

0.2mL/min (11MPa), 37°C, Single Quad. MS (ESI)

10µL (10µg/mL, 1% formic acid)