

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18


Cadenza CW-C18

Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Sulfasalazine (Antirheumatic)

スルファサラジン (抗リウマチ薬)

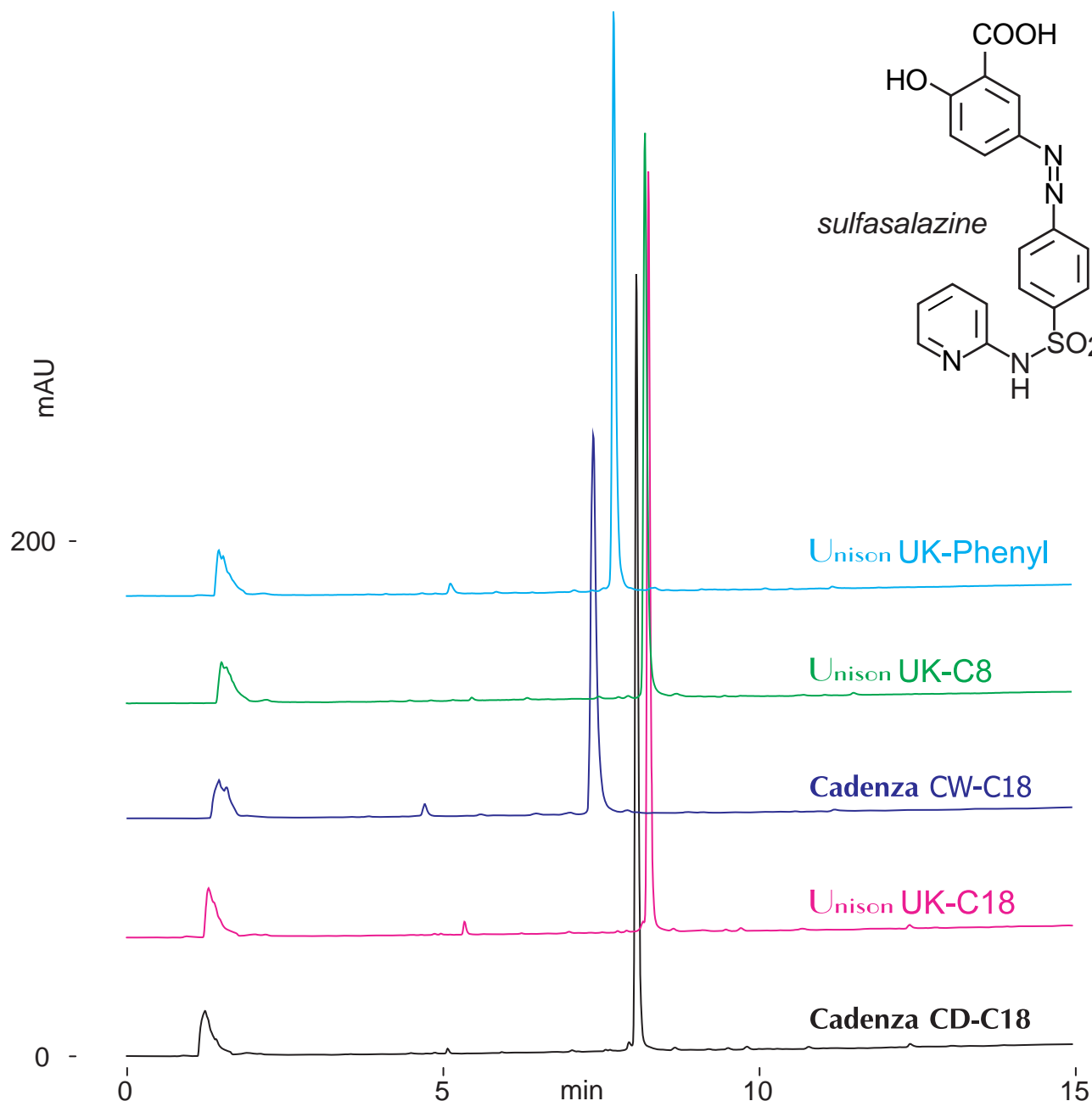

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

5-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.1mg/mL, dimethyl sulfoxide / formic acid = 100 / 1)