

Unison UK-Phenyl Cadenza CW-C18


Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Sulfasalazine (Antirheumatic)

スルファサラジン (抗リウマチ薬)

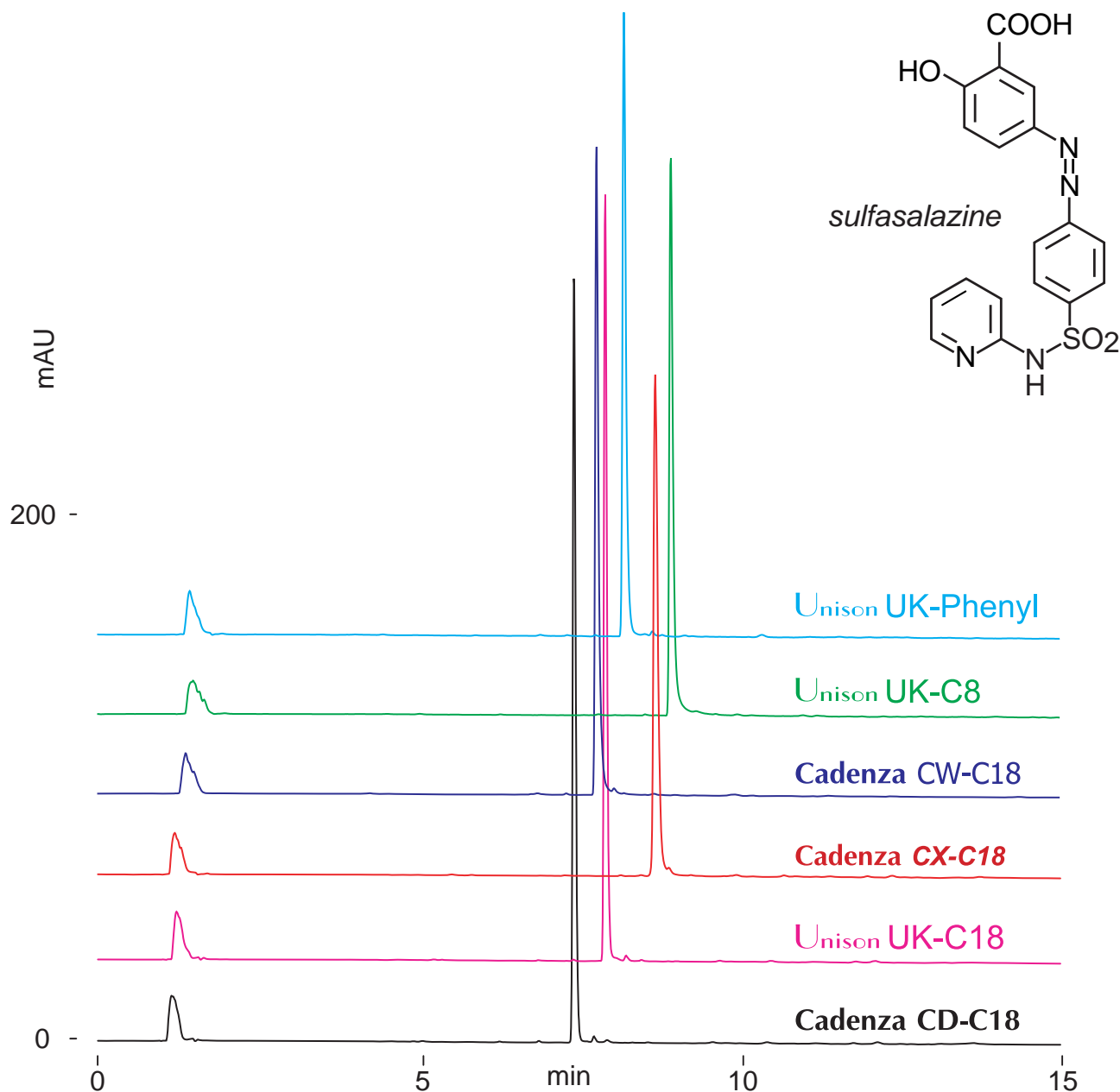

100 x 4.6 mm

A: acetonitrile / 100mM ammonium formate = 5 / 95

B: acetonitrile / 100mM ammonium formate = 90 / 10

0-65%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.1mg/mL, 100mM ammonium formate / dimethyl sulfoxide = 10 / 90)