



Scherzo SW-C18
 Scherzo SM-C18
 Scherzo SS-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Dobesilate (Vasoprotectives)

ドベシル酸 (血管保護薬)

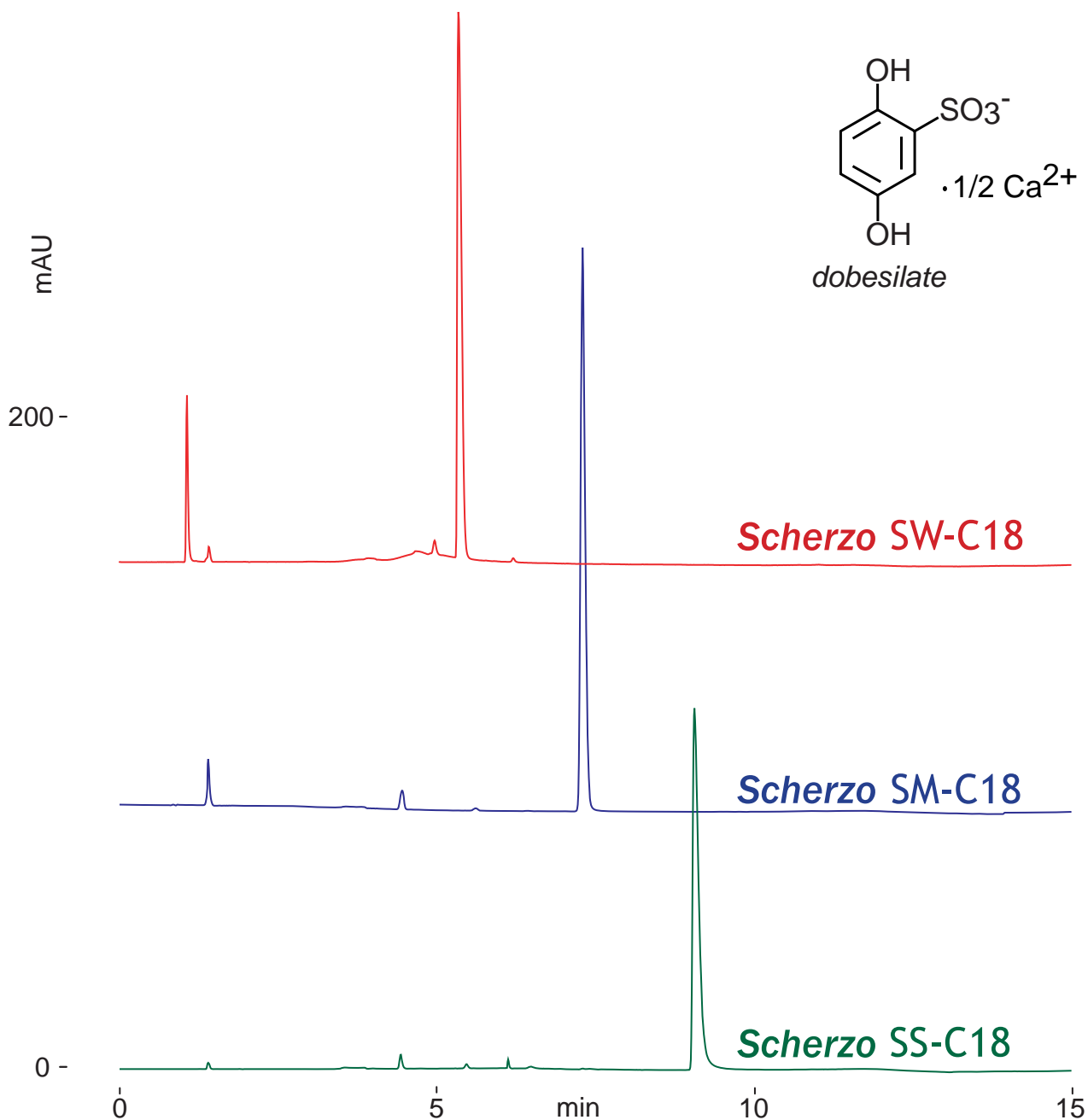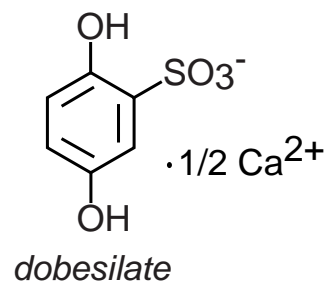

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: 100mM ammonium formate / acetonitrile = 50 / 50

0%B (0-1min), 0-100%B(1-10min), 100%B(10-15min)

1.0mL/min (9MPa), 37°C, 300nm

10µL (0.2mg/mL, water / formic acid = 100 / 1)