


Unison UK-C18
Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Melatonin related compounds

メラトニン関連化合物

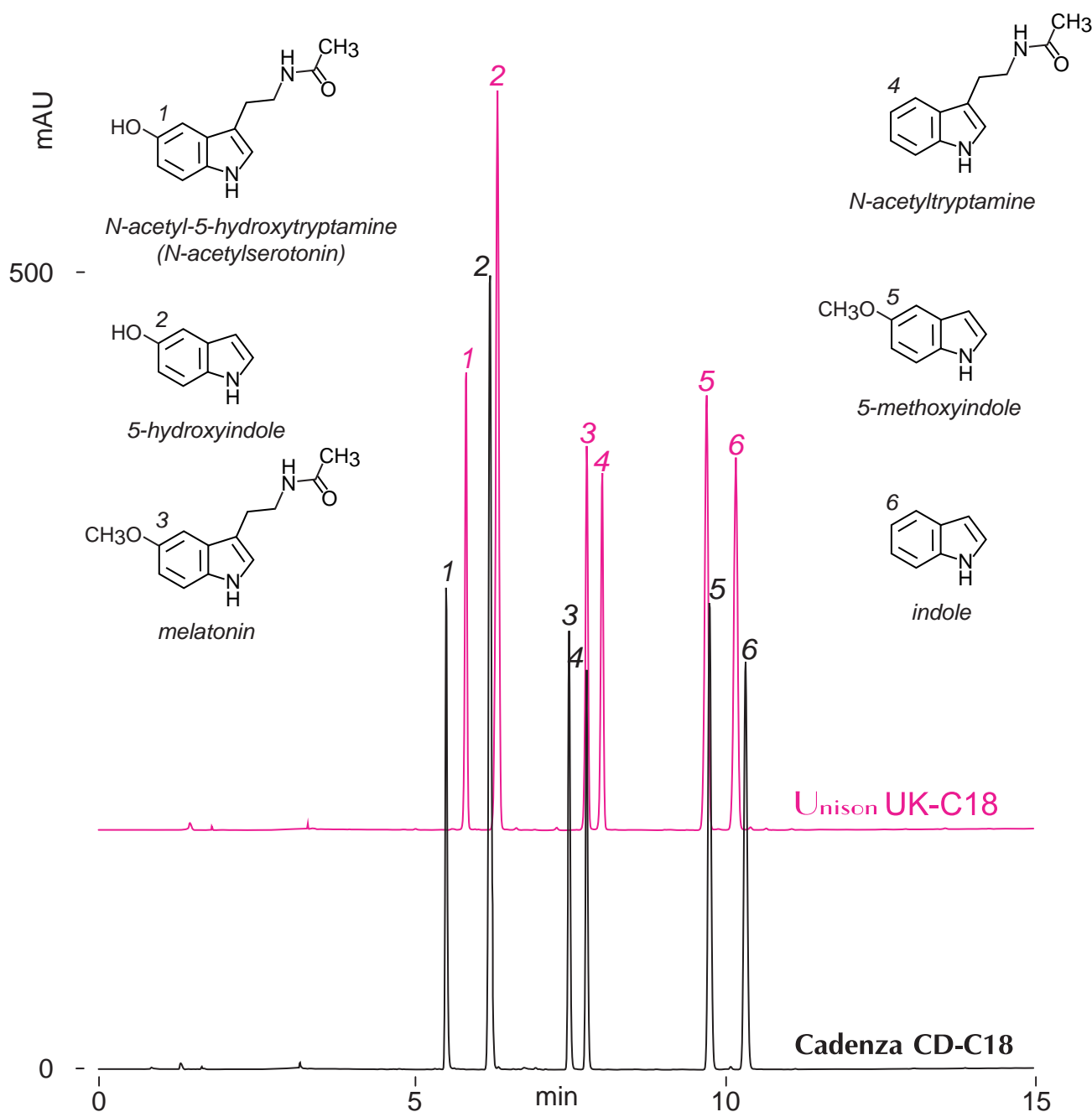

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.15mg/mL, water / acetonitrile / formic acid = 70 / 30 / 1)