

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Bezafibrate (Hypolipidemic agent)

ベザフィブラート (高脂血症治療薬)

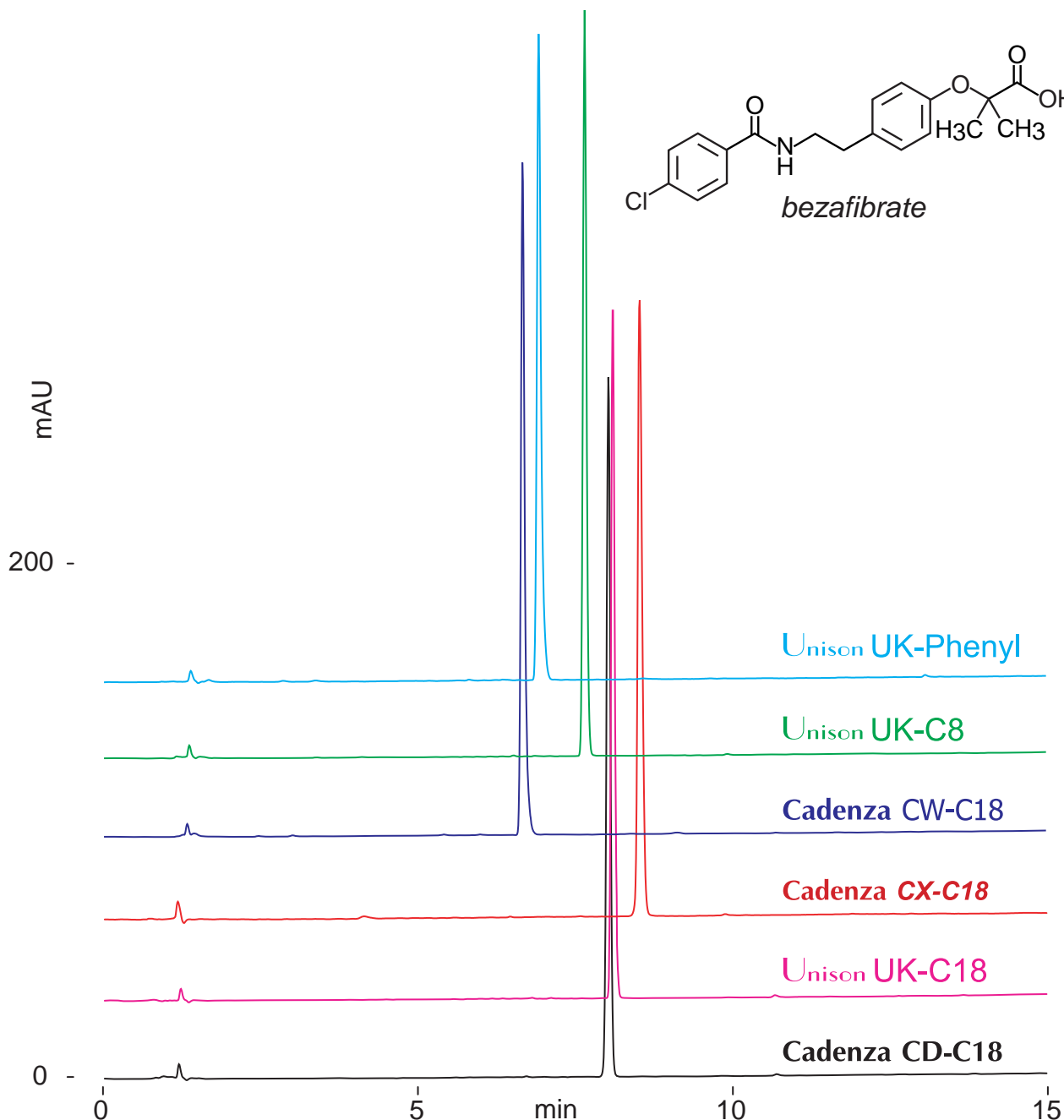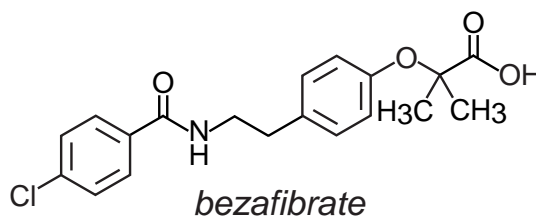

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

25-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.25mg/mL, , acetonitrile / water / formic acid = 50 / 50 / 1)