

Unison UK-Phenyl Cadenza CW-C18



Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Trimebutine (Stomach medicine)

トリメブチン (胃腸薬)

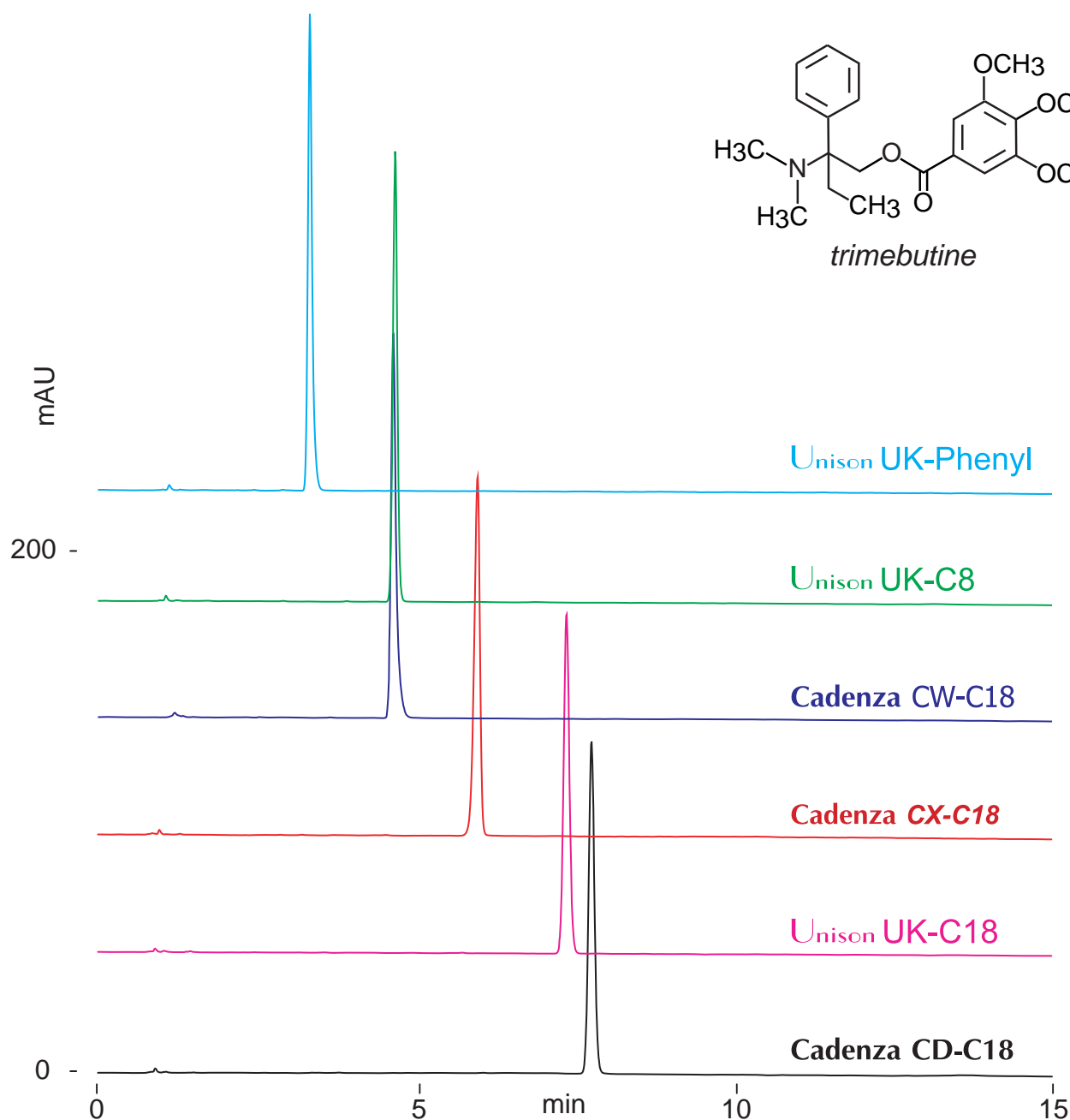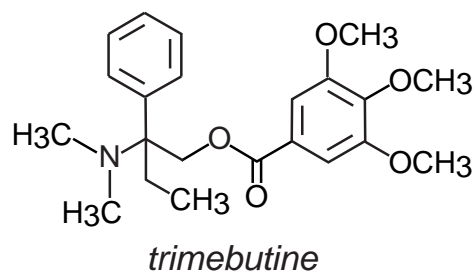

100 x 4.6 mm

A: acetonitrile / 100mM ammonium formate = 50 / 50

B: acetonitrile / 100mM ammonium formate = 90 / 10

10-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.15mg/mL, acetonitrile / 100mM ammonium formate = 50 / 50)