


Presto FF-C18

250 x 4.6 mm

HPLC COLUMN APPLICATION

Transferrin  
トランスフェリン

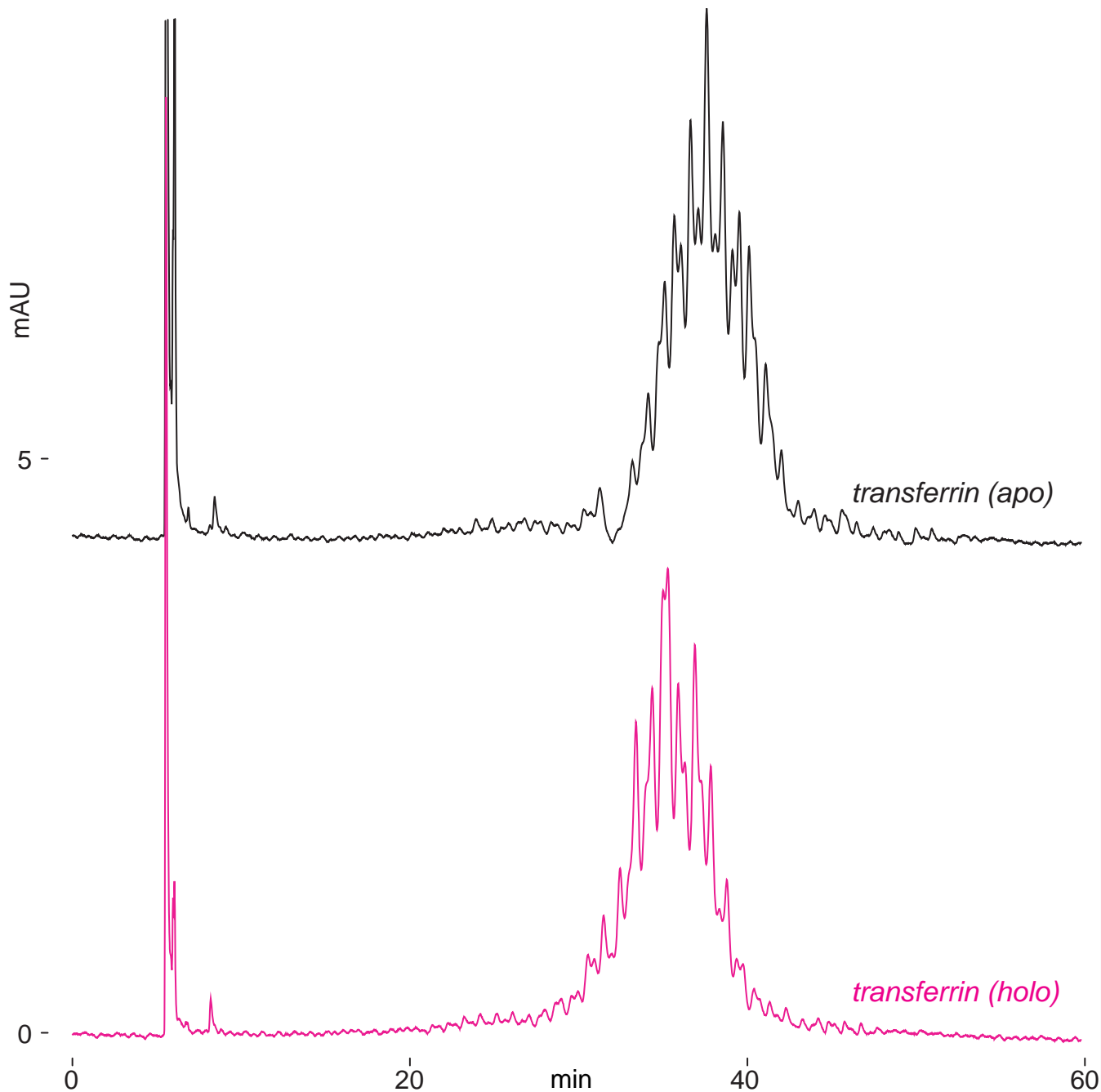

Presto FF-C18, 250 x 4.6 mm

A: water / acetonitrile / trifluoroacetic acid = 90 / 10 / 0.1

B: water / acetonitrile / trifluoroacetic acid = 10 / 90 / 0.1

34-38%B (0-60min)

0.3mL/min (22MPa), 37°C, 280nm

10µL (1mg/mL, water / trifluoroacetic acid = 100 / 0.1)