

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18



Cadenza CW-C18

Cadenza CD-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

## Carboxyphenylboronic acid


カルボキシフェニルボロン酸

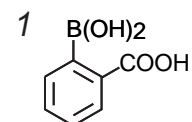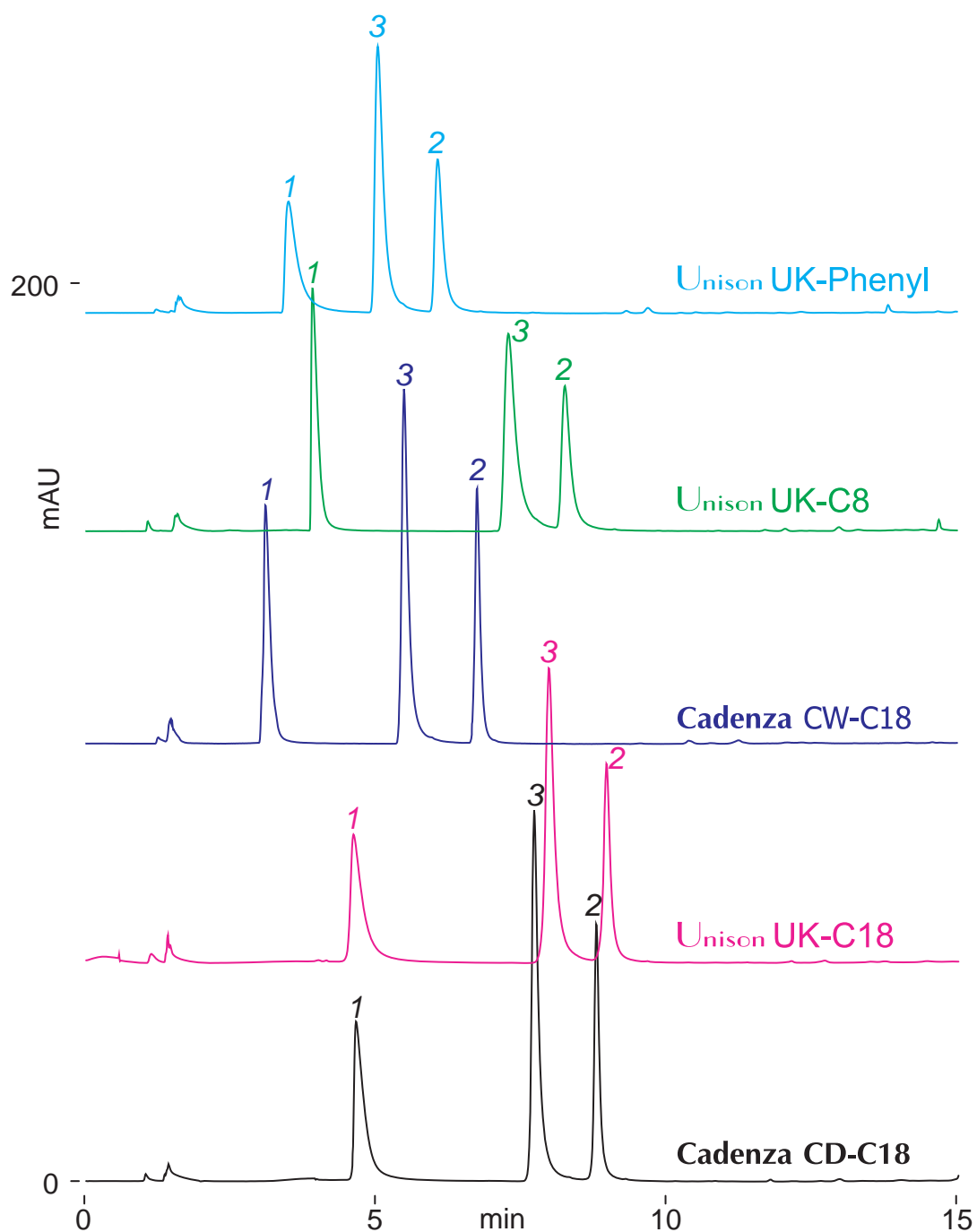

2-carboxyphenyl boronic acid

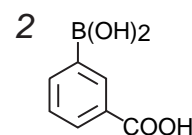

3-carboxyphenyl boronic acid

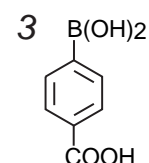

4-carboxyphenyl boronic acid

100 x 4.6 mm

A: water / formic acid = 100 / 0.5

B: water / acetonitrile / formic acid = 10 / 90 / 0.5

0-40%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.5mg/mL, water / DMSO / formic acid = 50 / 50 / 1)