

Unison UK-Phenyl Cadenza CX-C18


Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Bosentan, Macitentan (pulmonary hypertension drug)

ボセンタン, マシテンタン (肺高血圧治療薬)

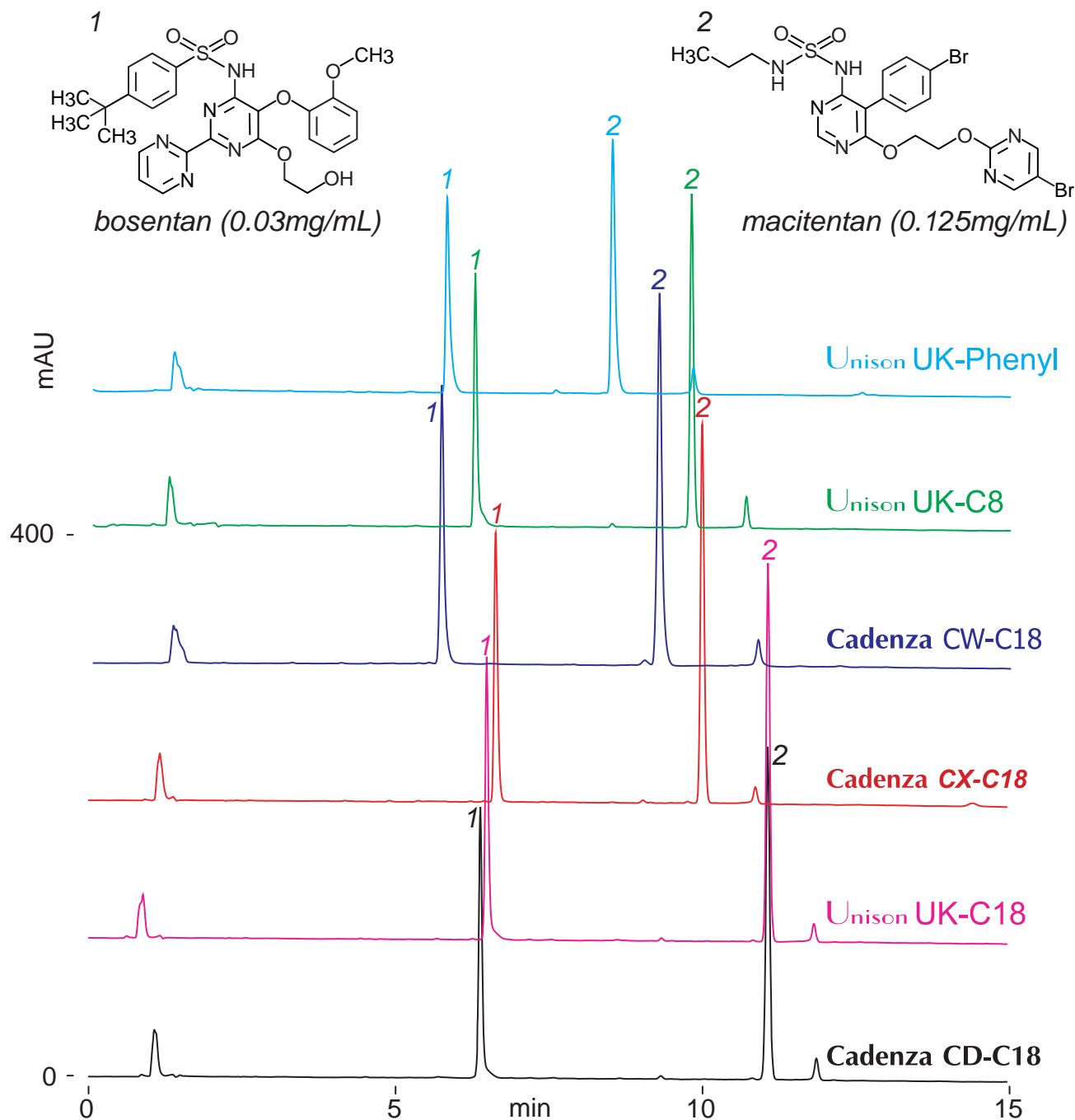

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

25-100%B (0-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.03-0.125mg/mL, DMSO / 100mM ammonium formate = 50 / 50)