



Scherzo SW-C18  
 Scherzo SM-C18  
 Scherzo SS-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Folinic acid (Antifolate)

ホリン酸 (葉酸代謝拮抗薬)

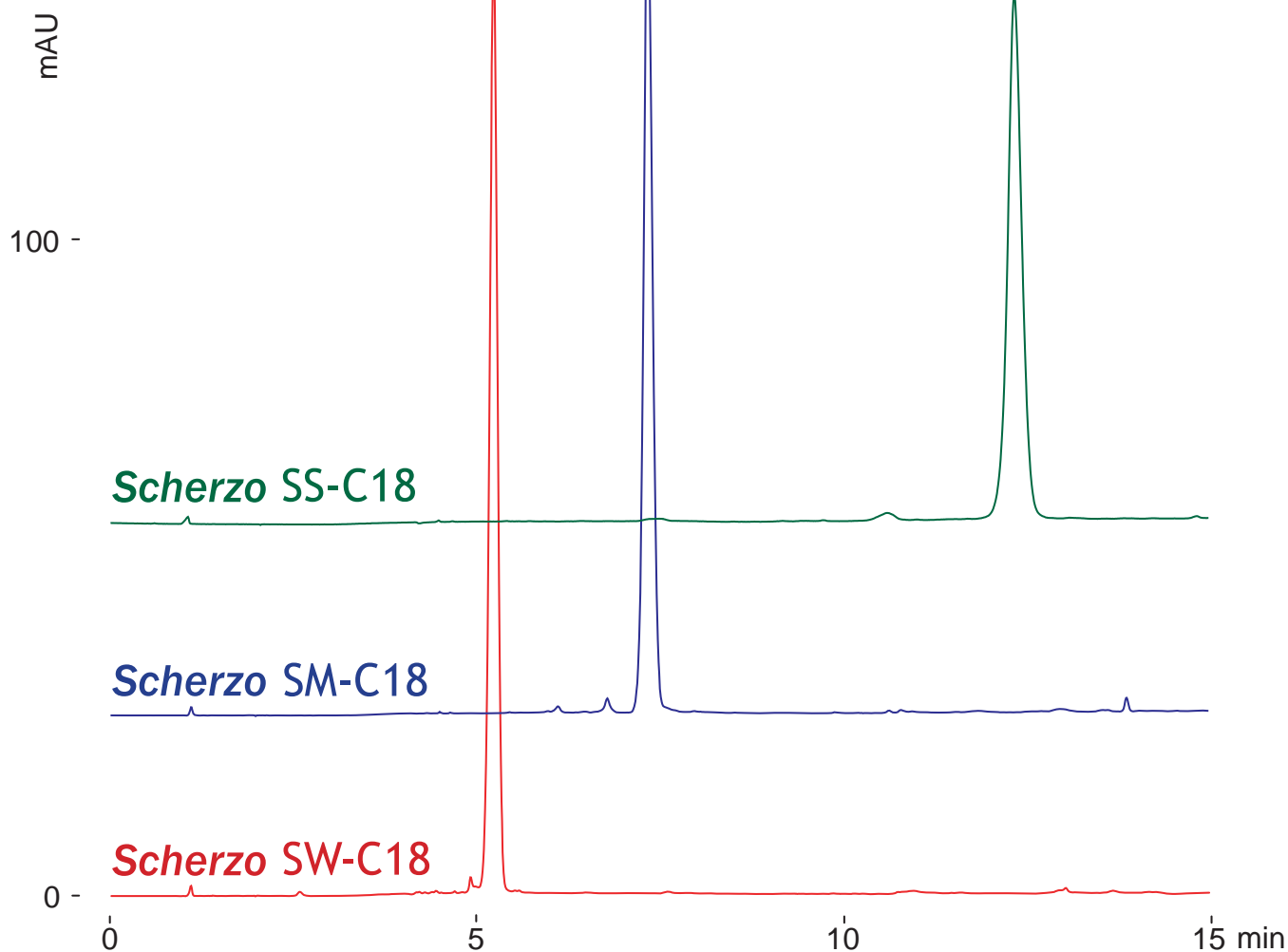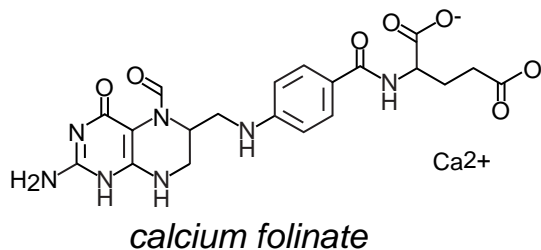

100 x 4.6 mm  
 A: 10mM ammonium formate  
 B: 100mM ammonium formate / acetonitrile = 50 / 50  
 0%B (0-1min), 0-100%B (1-15min)  
 1.0mL/min (8MPa), 37°C, 260nm  
 10µL (0.15mg/mL, 100mM ammonium formate )