


Cadenza CX-C18

150 x 3 mm

HPLC COLUMN APPLICATION

Loratadine (Antiallergic agent)

ロラタジン (抗アレルギー薬)

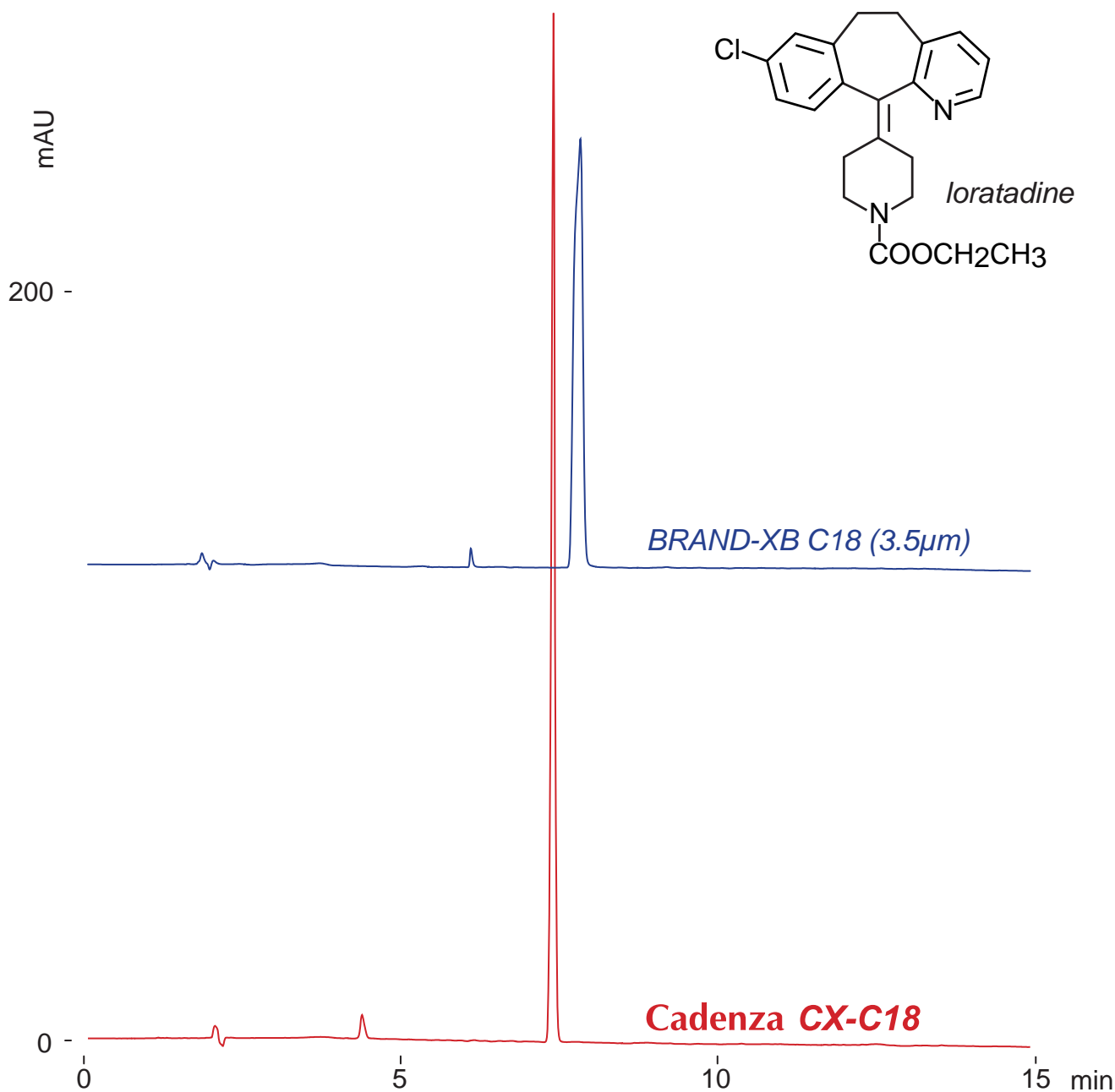

150 x 3 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-100%B (0-15min)

0.4mL/min (10MPa), 37°C, 260nm

2µL (0.5mg/mL, water / acetonitrile / formic acid = 50 / 50 / 1)