
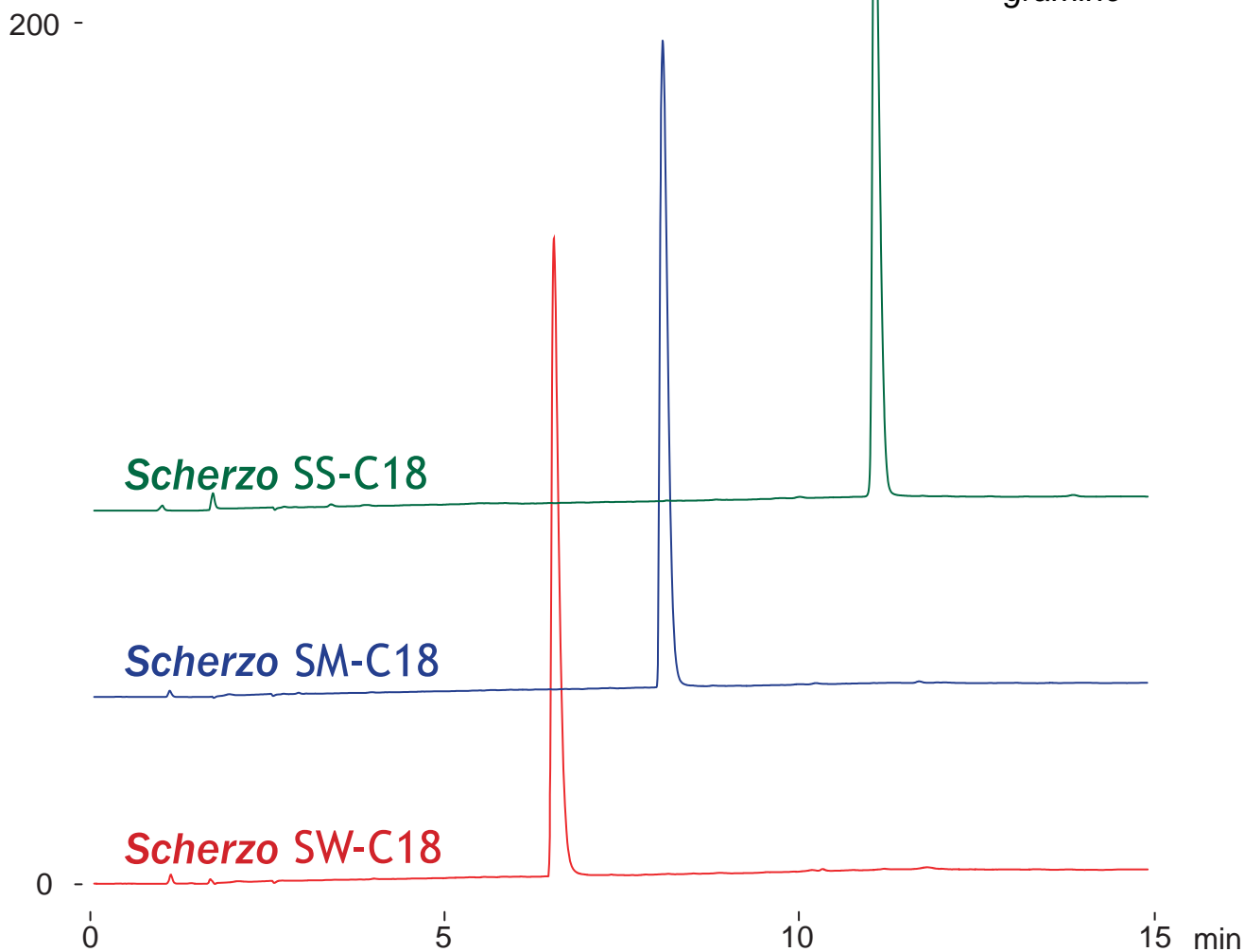


Scherzo SW-C18  
 Scherzo SM-C18  
 Scherzo SS-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

**Gramine (Indole alkaloid)**  
 グラミン (インドールアルカロイド)

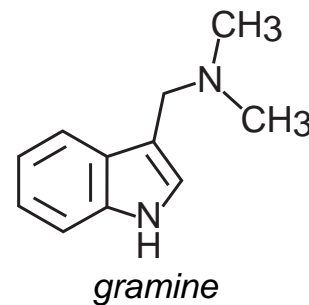

100 x 4.6 mm

A: 10mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 50 / 50

0-100%B (0-10min), 100%B(10-15min)

1.0 mL/min (12MPa), 37°C, 260nm

10µL (0.12mg/mL, 100mM ammonium formate / acetonitrile = 67 / 33 )