

Unison UK-Phenyl Cadenza CX-C18


Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

Isoflavones (Phytoestrogens)


イソフラボン類 (植物エストロゲン)

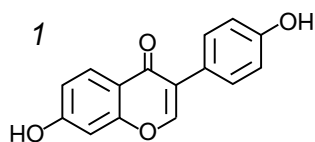


daidzein
(0.048mg/mL)

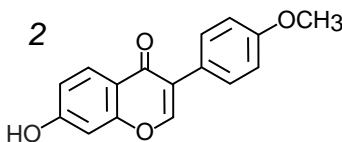

formononetin
(0.029mg/mL)

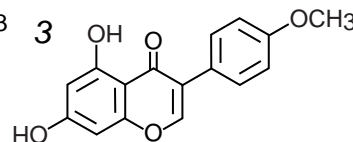

biochanin A
(0.045mg/mL)

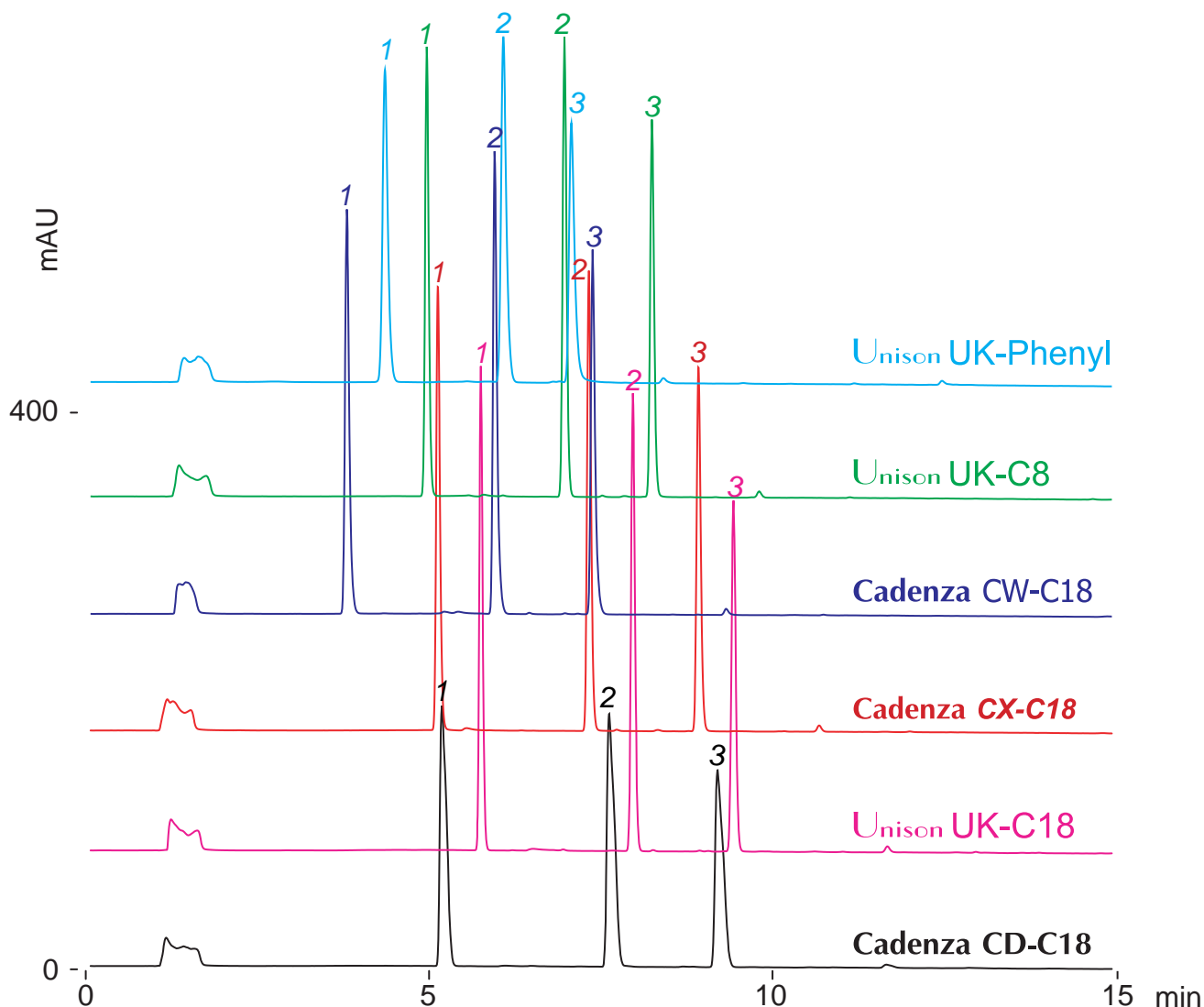

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.029-0.045mg/mL, 100mM ammonium formate / DMF = 10 / 90)