



Unison UK-Phenyl Cadenza CX-C18

Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Brucine (Indole alkaloid)
 ブルシン (インドールアルカロイド)

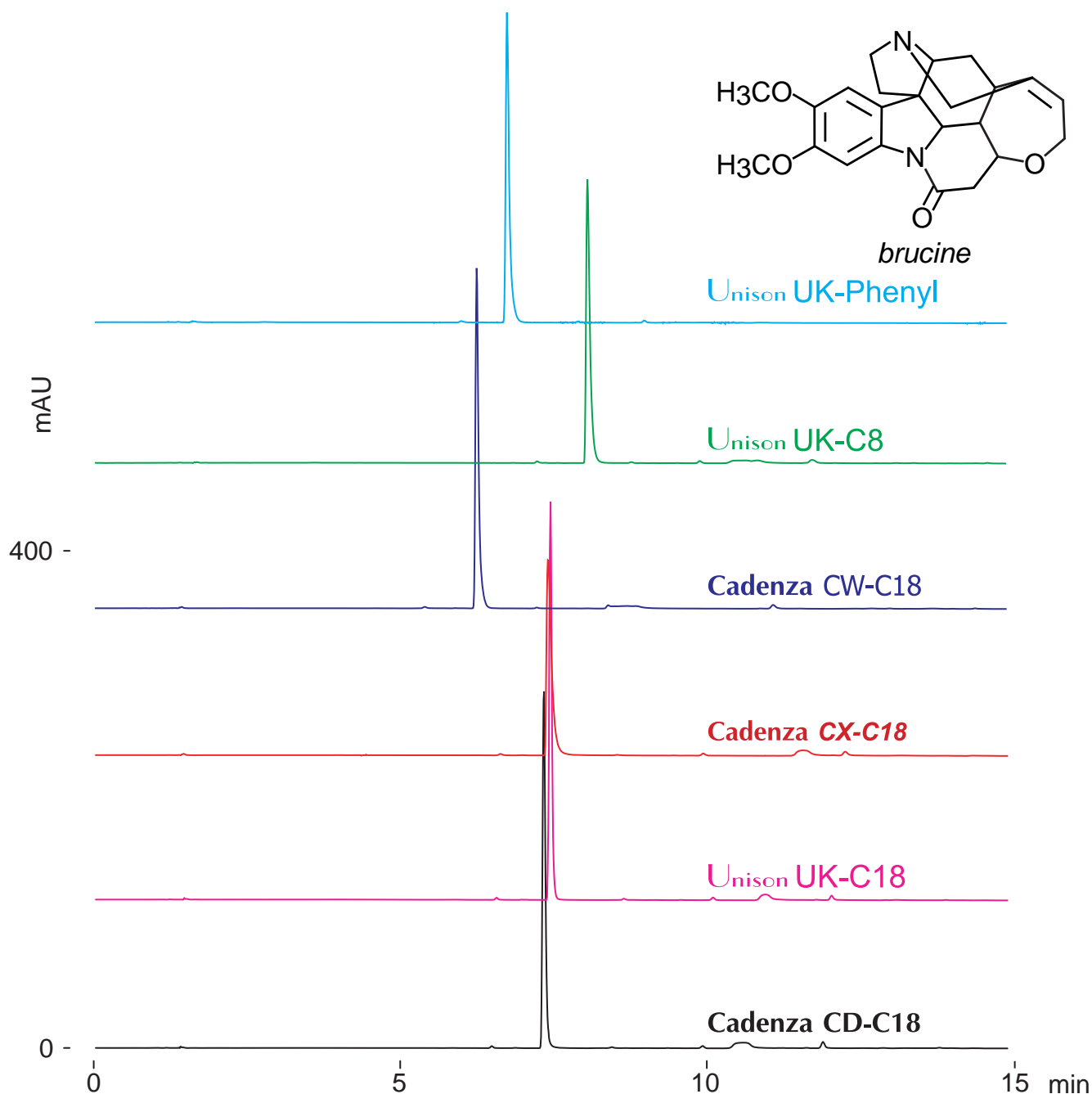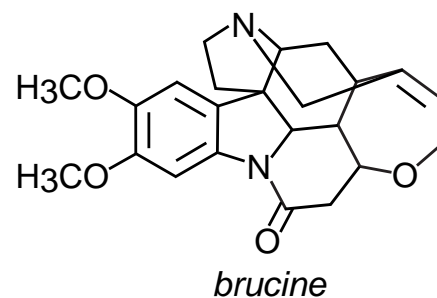

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

0%B (0-1min), 0-50%B (1-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.1mg/mL, 100mM ammonium formate / acetonitrile = 75 / 25)