

Cadenza CX-C18

Intrada WP-RP


Presto FF-C18

150 x 4.6 mm


HPLC COLUMN APPLICATION

## Polyhexamethylene biguanide, Benzalkonium


ポリヘキサメチレンビグアナイド、ベンザルコニウム

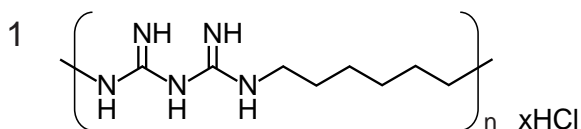

polyhexamethylene biguanide hydrochloride  
(PHMB, 1000ppm)

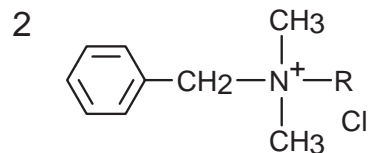

benzalkonium chloride  
(500ppm)

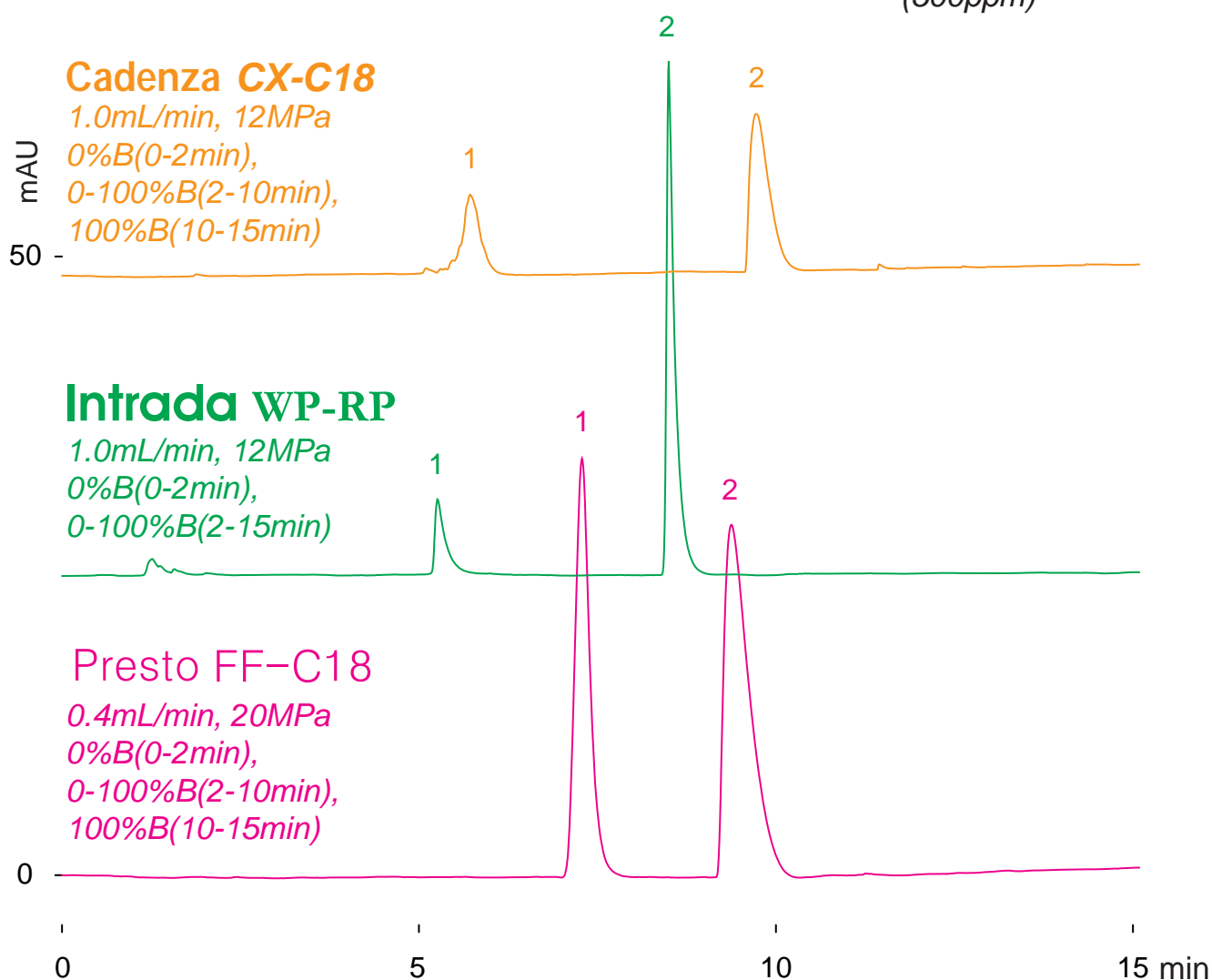

150 x 4.6 mm

A: water / acetic acid = 100 / 1

B: acetonitrile / acetic acid = 100 / 1

37°C, ELSD (spray chamber 70°C, drift chamber 100°C)

10µL (500-1000ppm, water)