

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Amygdalin (Cyanogenic glycoside)

アミグダリン (青酸配糖体)

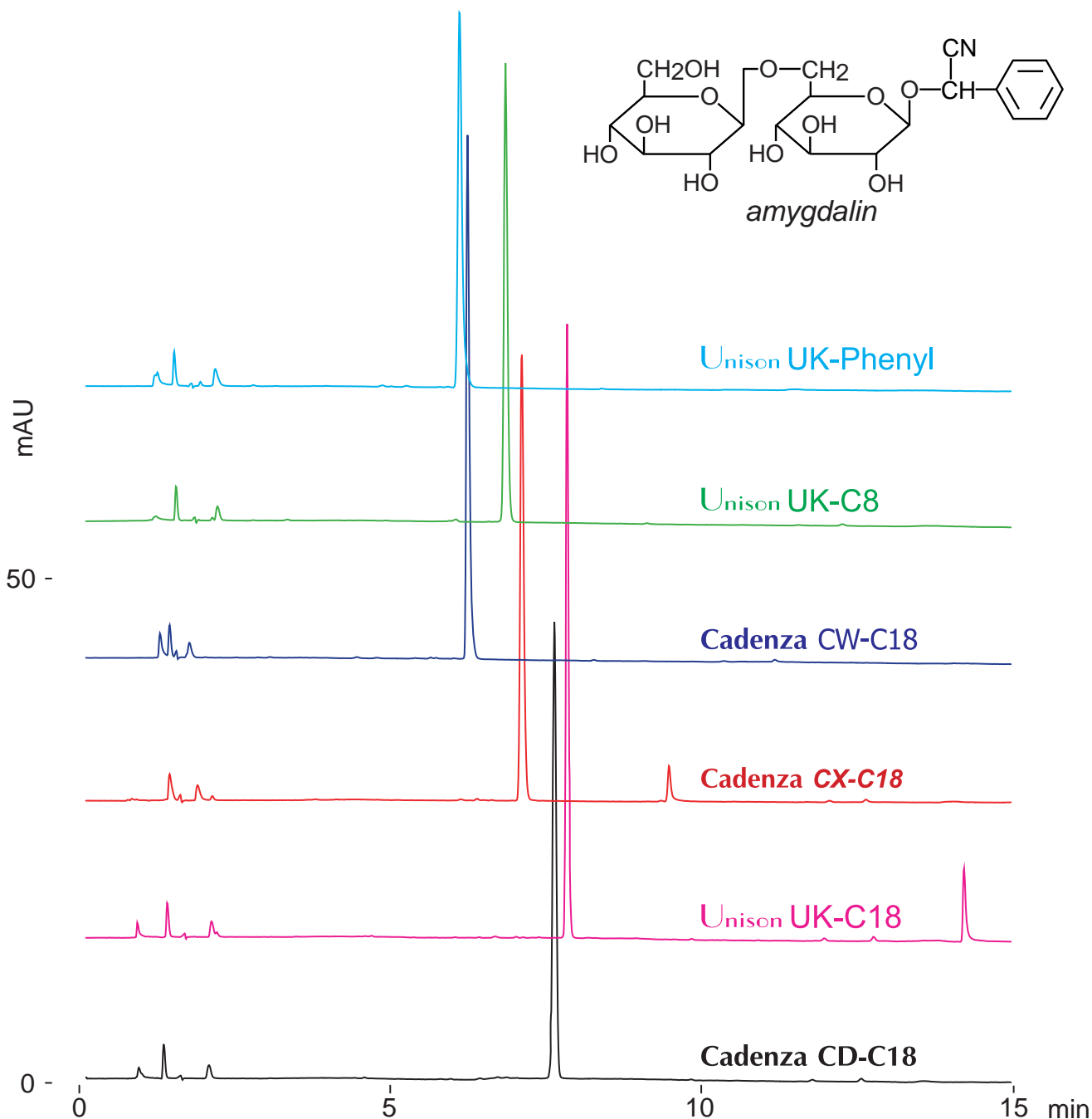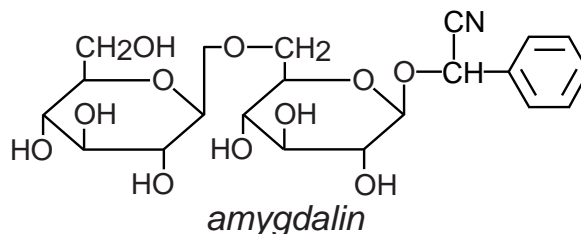

100 x 4.6 mm

A: water / acetonitrile / formic acid = 95 / 5 / 0.1

B: water / acetonitrile / formic acid = 50 / 50 / 0.1

0%B (0-1min), 0-60%B (1-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (1mg/mL, water / formic acid = 100 / 1)