




Scherzo SM-C18

250 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Neomycin B and C

ネオマイシン B及び Cの LC-MS分析

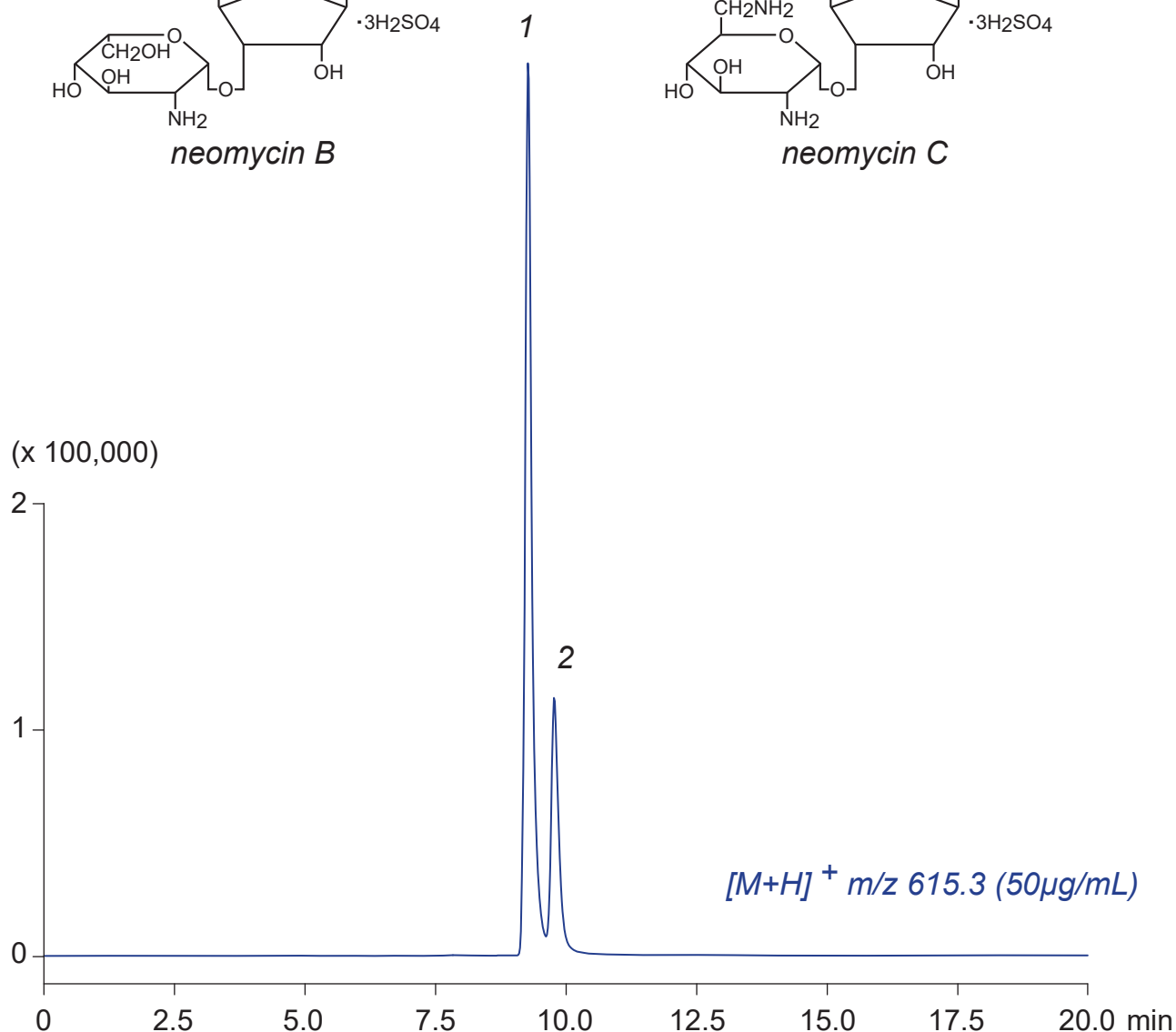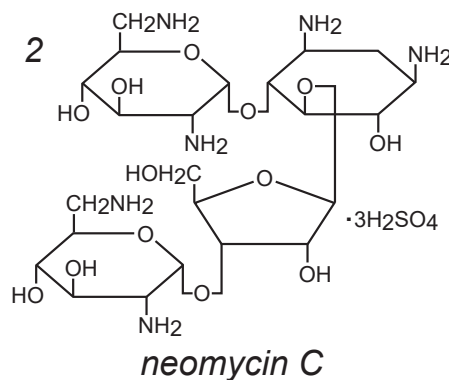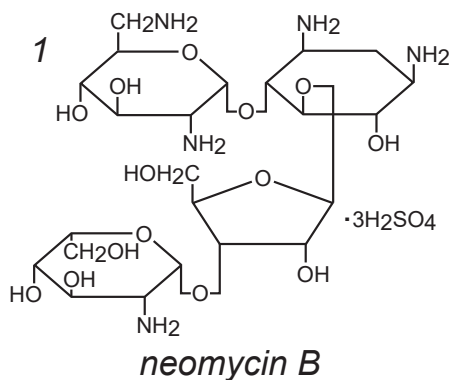

Scherzo SM-C18, 250 x 3 mm

A: 10mM ammonium formate

B: acetonitrile / water / formic acid = 50 / 50 / 0.1

5-40%B (0-20min)

0.4mL/min (14MPa), 37°C, Single Quad. MS (ESI)

10µL (50µg/mL, 100mM ammonium formate)