


Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Methylcobalamin (Coenzyme form of Vitamin B12)
 メチルコバラミン (補酵素型ビタミン12)

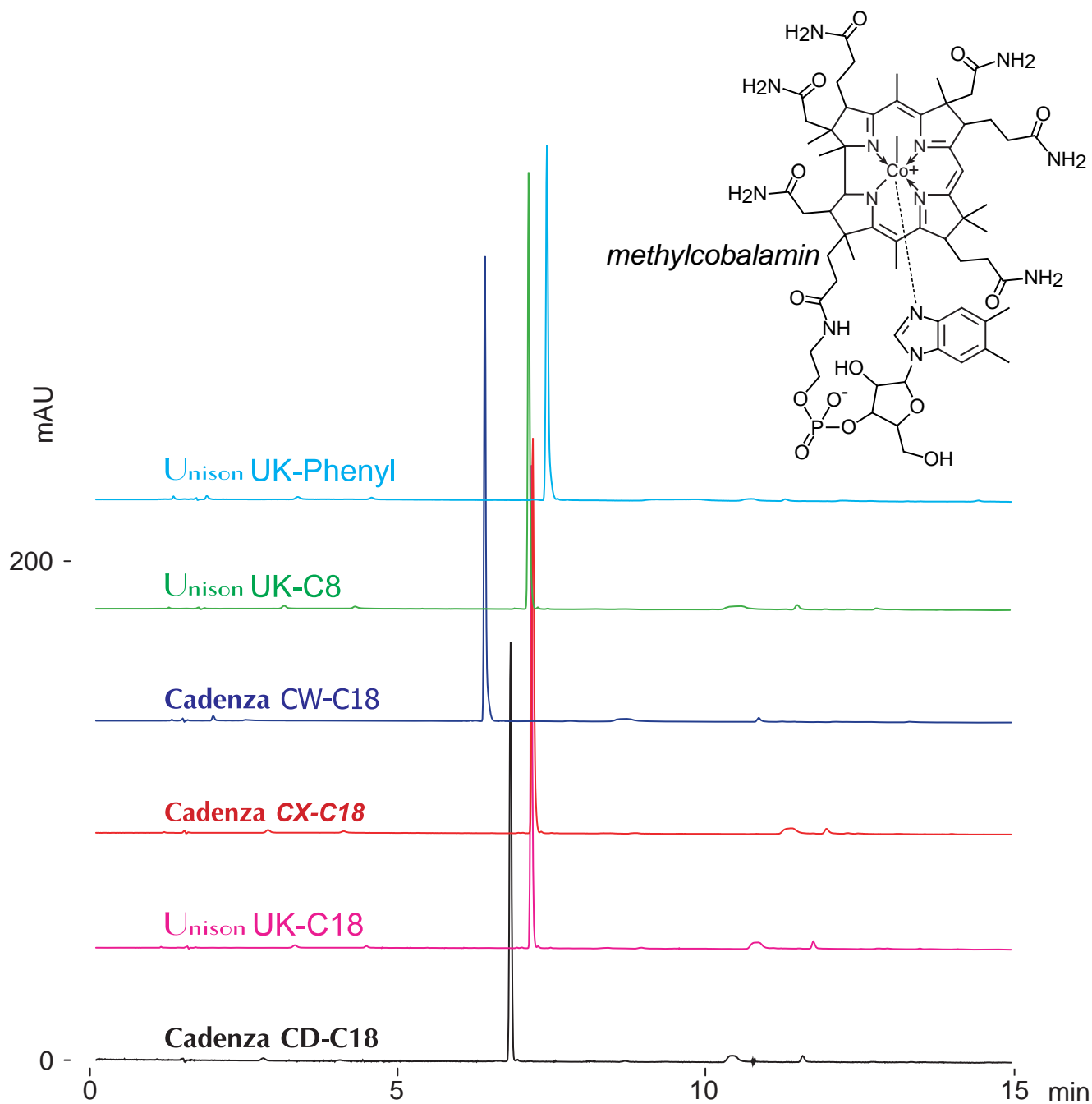

100 x 4.6 mm
 A: 100mM ammonium formate / acetonitrile = 90 / 10
 B: 100mM ammonium formate / acetonitrile = 10 / 90
 0%B (0-1min), 0-50%B (1-15min)
 1.0mL/min (9MPa), 37°C, 260nm
 10µL (0.15mg/mL, 100mM ammonium formate)