


**Intrada Organic Acid**

150 x 2 mm

**HPLC COLUMN APPLICATION**

**LC-MS: Dioctyl sulfosuccinate (Duble-tail anionic surfactant)**

ジオクチルスルホコハク酸 (二鎖型アニオン界面活性剤) の LC-MS 分析

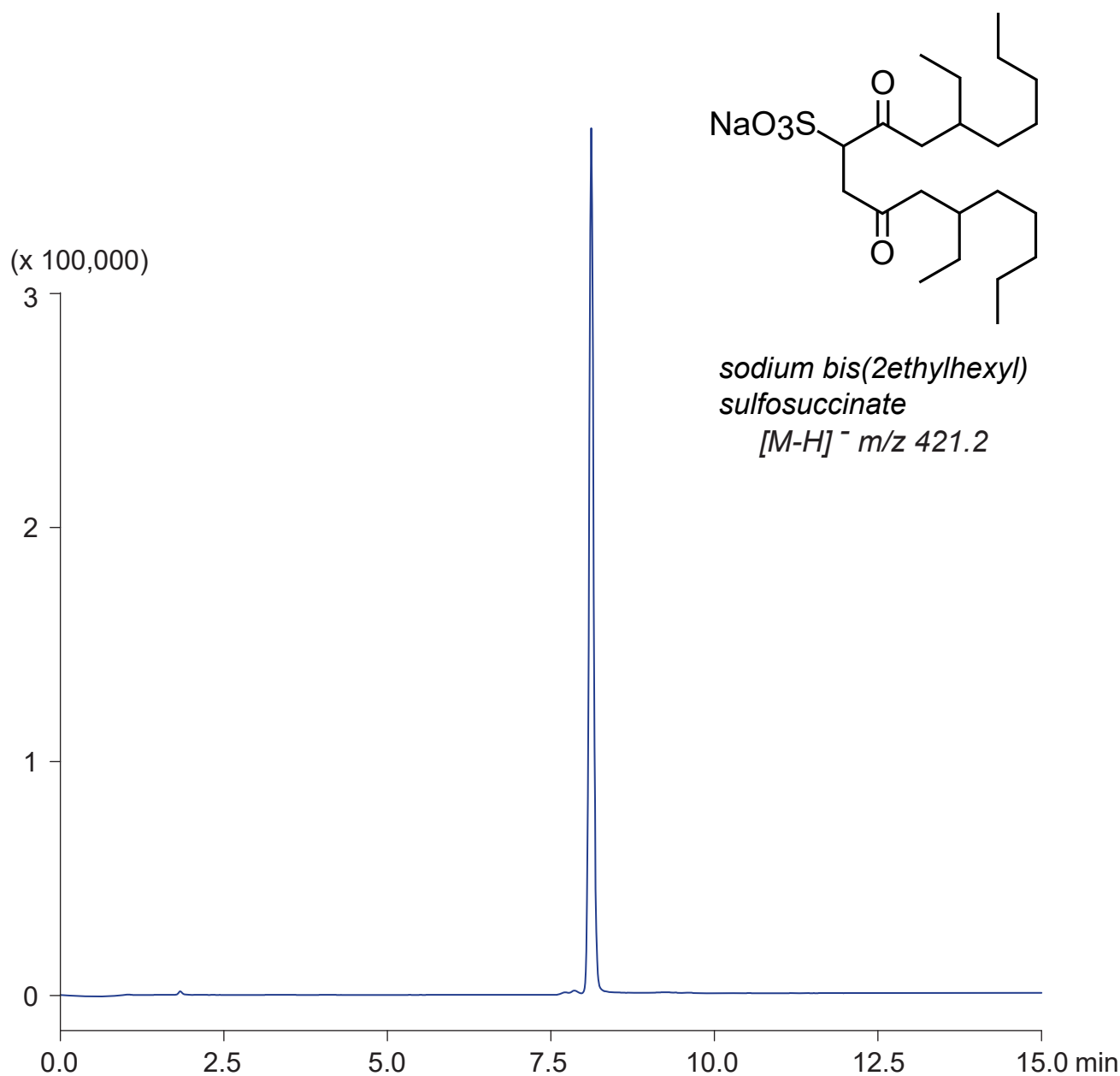

*Intrada Organic Acid, 150 x 2 mm*

*A: acetonitrile / water / formic acid = 30/ 70 / 0.1*

*B: acetonitrile / 100mM ammonium formate = 60/ 40*

*0-100%B (0-10min), 100%B (10-15min)*

*0.2mL/min (13MPa), 37°C, Single Quad. MS (ESI)*

*5uL (50ng/mL, ACN / water / HCOOH = 50 / 50 / 1)*