


**Scherzo SM-C18 MF**

100 x 2 mm



**HPLC COLUMN APPLICATION**

**LC-MS/MS: Glyphosate, N-acetylglyphosate**

グリホサート, N-アセチルグリホサートの LC-MS/MS 分析

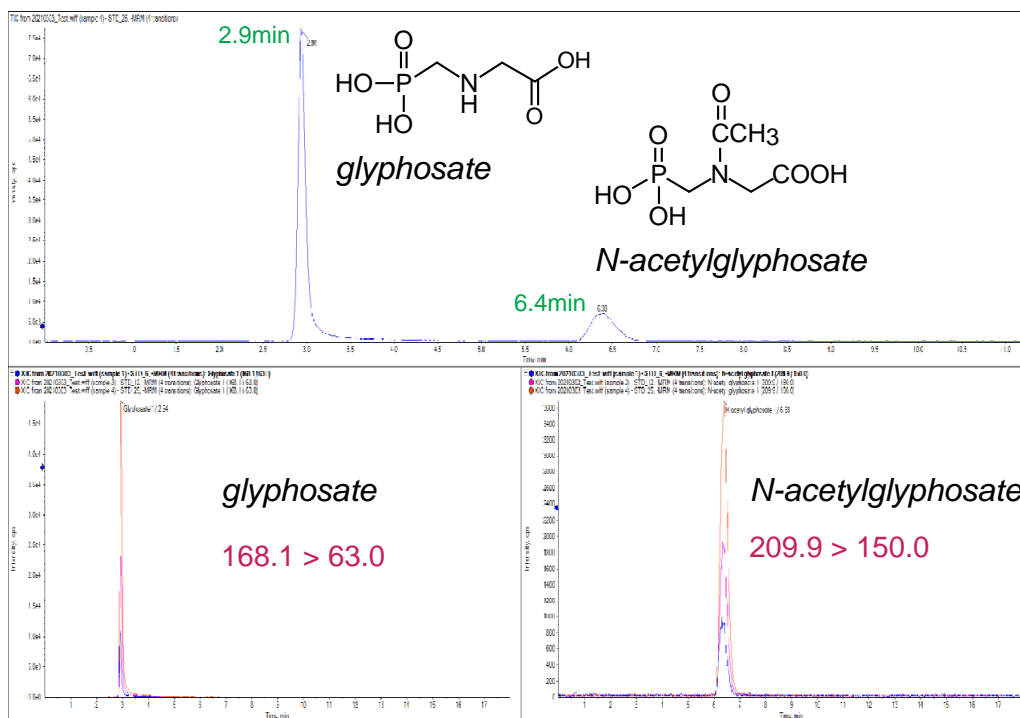

**Glyphosate STD, N-acetyl glyphosate 6, 12, 25 ppb cal curve**

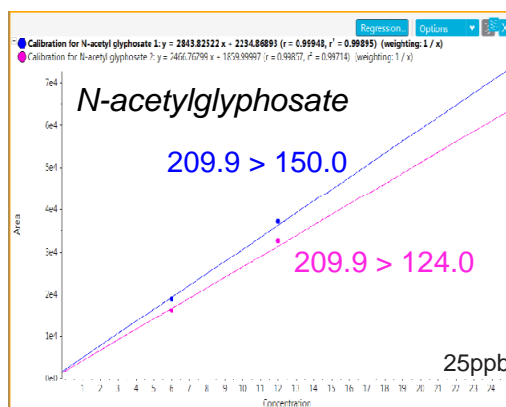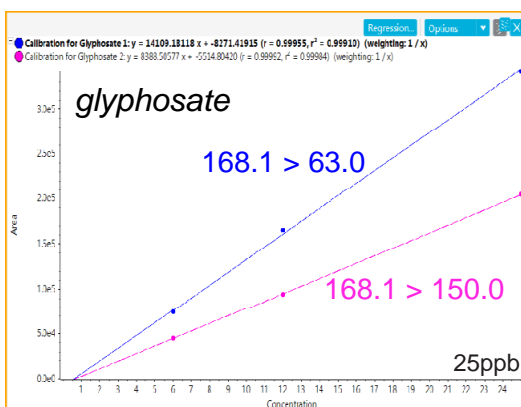

Scherzo SM-C18 MF, 100 x 2 mm

A: 0.1% formic acid in 20% methanol

B: 1% formic acid in 20% methanol

0%B (0-1min), 0-100%B (1-4min), 100%B (4-11min),

100-0%B (11-11.1min), 0%B (11.1-18min)

0.25mL/min, 40°C, Triple Quad. MS (ESI, Negative)

5µL (6, 12, 25ppb)

Courtesy of PEAKMAN SP, KOREA